

ST™

by MEGA'GEN

11° internal hex connection

KnifeThread combined
with double thread

S-L-A surface
treatment for safety &
better bone growth





Key Advantages

11° internal hex connection

KnifeThread combined with double thread

S-L-A surface treatment for safety & better bone growth



04 Characteristics & Advantages

04 I. Features

05 Fixture Product

05 I. Dimension

06 II. Size

07 Cover Screw & Healing Abutment

10 Abutment & Prosthetic Options

10 I. Fixture Level Prosthesis

17 1. Fixture Level Prosthesis_Digital

27 II. Abutment Level Prosthesis

27 1. Solid Abutment & Components

29 2. Octa Abutment & Components

35 3. Multi-unit Abutment & Components

47 ST Kit

47 I. Surgical Kit

52 Clinical Case


Characteristics & Advantages

I. Features

AnyRidge® technology maintains
100% compatibility with existing systems


More economical

11° internal connection with high compatibility

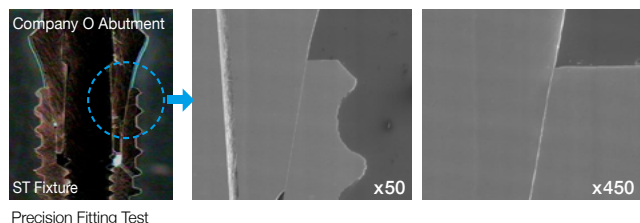



KnifeThread® is added to Double Thread

Thread design is very important for surgery as well as threading. The KnifeThread® of ST™ system has powerful self-threading capability that ensures initial stability, heightens resistance to compressive force, and minimizes the occurrence of shear force, leading to excellent load distribution.

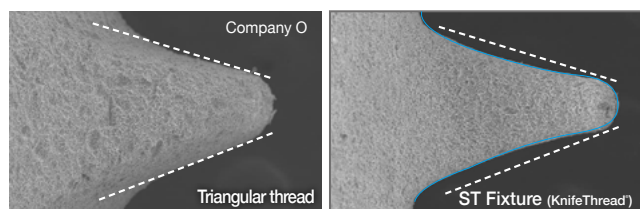

S-L-A surface treatment based on advanced technology

- S-L-A surface morphology combines macro & micro roughness via blasting & etching processes
- Proven to be safer, as sharp edges can crack & become foreign objects during implant placement
- Superior bone formation capability

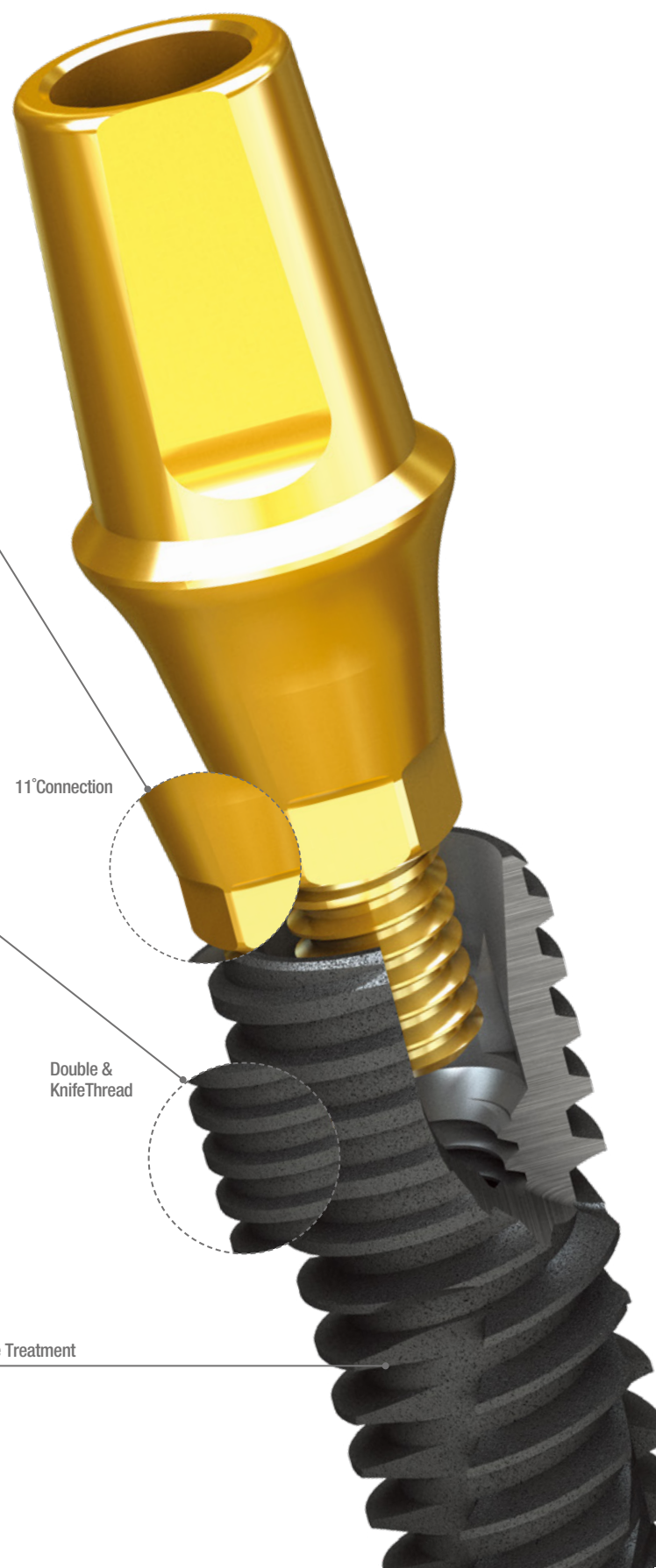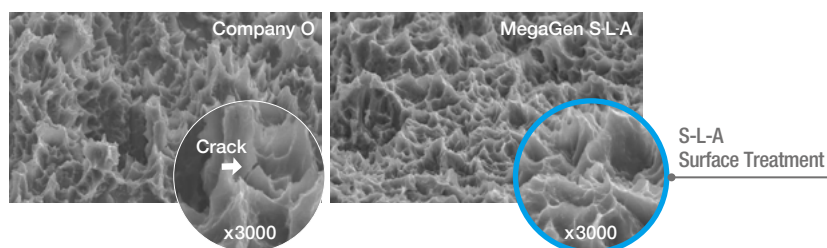


Fixture Product

I. Dimensions


Platform Diameter & Female Screw

Type	Diameter	Hex(mm)	Thread
Mini	Ø3.5	2.1	M1.6 (M2.0)
Regular	Ø4.0	2.5	M2.0
	Ø4.5		
	Ø5.0		

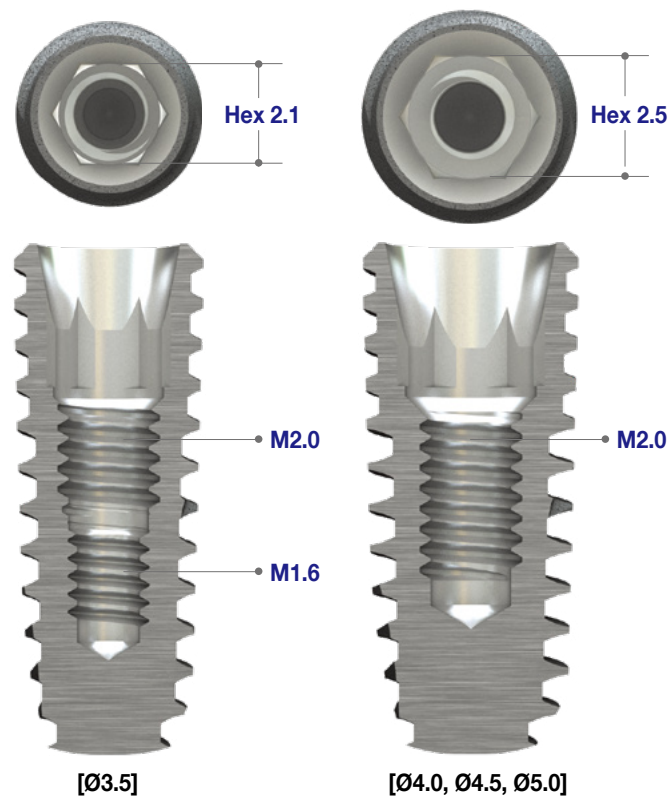

Fixture Size Variation

Diameter	Widest thread Diameter	Length(mm)	Platform Diameter
Ø3.5	Ø3.7	8.5 / 10.0 / 11.5 / 13.0 / 15.0	Ø3.1
Ø4.0	Ø4.2	7.0 / 8.5 / 10.0 / 11.5 / 13.0 / 15.0	Ø3.6
Ø4.5	Ø4.6	7.0 / 8.5 / 10.0 / 11.5 / 13.0 / 15.0	Ø4.0
Ø5.0	Ø5.1	7.0 / 8.5 / 10.0 / 11.5 / 13.0 / 15.0	Ø4.5

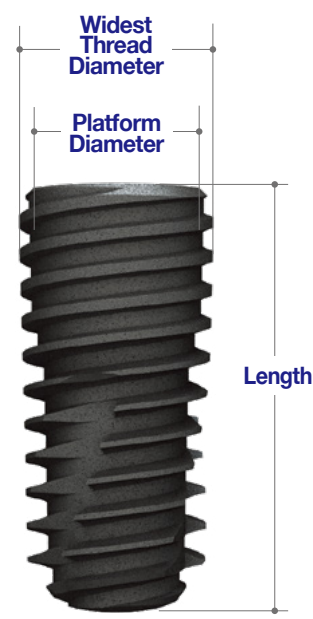

II. Sizes

Mini Ø3.5

- Cover Screw (STCSS) included

Diameter	Length(mm)	Ref.C
Ø3.5	8.5	ST3508C
	10.0	ST3510C
	11.5	ST3511C
	13.0	ST3513C
	15.0	ST3515C



Regular Ø4.0

- Cover Screw (STCSR) included

Diameter	Length(mm)	Ref.C
Ø4.0	7.0	ST4007C
	8.5	ST4008C
	10.0	ST4010C
	11.5	ST4011C
	13.0	ST4013C
	15.0	ST4015C



Regular Ø4.5

- Cover Screw (STCSR) included

Diameter	Length(mm)	Ref.C
Ø4.5	7.0	ST4507C
	8.5	ST4508C
	10.0	ST4510C
	11.5	ST4511C
	13.0	ST4513C
	15.0	ST4515C



Regular Ø5.0

- Cover Screw (STCSR) included

Diameter	Length(mm)	Ref.C
Ø5.0	7.0	ST5007C
	8.5	ST5008C
	10.0	ST5010C
	11.5	ST5011C
	13.0	ST5013C
	15.0	ST5015C



Cover Screw & Healing Abutment

Cover Screw

* Included in fixture package

- Hand driver (1.2 Hex) is used
- Used for submerged-type surgery
- Internal structure of fixture
- Recommended torque: by hand (5 - 8Ncm)


Height (mm)	Type	Ref.C
0.4	Mini	STCSS
	Regular	STCSR



Healing Abutment

- Hand Driver (1.2 Hex) is used
- Used for submerged-type surgery or two-stage surgery
- Various diameters & heights for different cases
- Forms appropriate emergence profile during gingival healing
- Recommended torque: by hand (5 - 8Ncm)

Diameter	Type	Height(mm)	Ref.C
Ø4.0	Mini	3.0	STHA403M
		4.0	STHA404M
		5.0	STHA405M
		7.0	STHA407M
		9.0	STHA409M
Ø4.5	Mini	3.0	STHA453M
		4.0	STHA454M
		5.0	STHA455M
		7.0	STHA457M
		9.0	STHA459M
Ø4.0	Regular	3.0	STHA403R
		4.0	STHA404R
		5.0	STHA405R
		7.0	STHA407R
		9.0	STHA409R
Ø4.5	Regular	3.0	STHA453R
		4.0	STHA454R
		5.0	STHA455R
		7.0	STHA457R
		9.0	STHA459R
Ø5.0	Regular	3.0	STHA503R
		4.0	STHA504R
		5.0	STHA505R
		7.0	STHA507R
		9.0	STHA509R
Ø6.0	Regular	3.0	STHA603R
		4.0	STHA604R
		5.0	STHA605R
		7.0	STHA607R
		9.0	STHA609R
Ø7.0	Regular	3.0	STHA703R
		4.0	STHA704R
		5.0	STHA705R
		7.0	STHA707R
		9.0	STHA709R

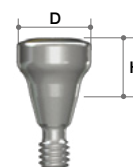

All in one package


(Anatomic type)

Height (mm)	Connection	Ref.C
4	Hex	C-AHA004HP
	Non-Hex	C-AHA004NP
5	Hex	C-AHA005HP
	Non-Hex	C-AHA005NP
7	Hex	C-AHA007HP
	Non-Hex	C-AHA007NP

Healing Abutment

(Anatomic type)

- Use with a Hand Driver(1.2 Hex).
- Abutment Screw included.H=4 AOHAS2004/ H=5 AOHAS2005/ H=7 AOHAS2007
- Used for non-submerged type surgery or for two stage surgery.
- Choose appropriate diameter and height of Healing Abutment according to situation.
- Helps to form suitable emergence profile during period of gingival healing.
- Recommend torque : by hand (5 - 8Ncm)

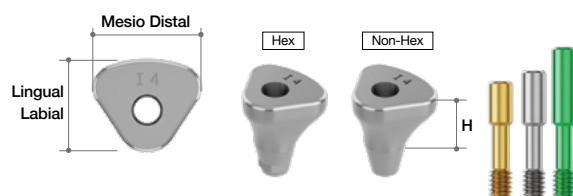

Type	MD (mm)	LL (mm)	Position No.	Height (mm)	Connection	Ref.C
Incisor	4.0	5.0	I1	4	Regular	AOHI40504T
				5		AOHI40505T
				7		AOHI40507T
	4.5	4.5	I2	4		AOHI45454T
				5		AOHI45455T
				7		AOHI45457T
	6.0	5.0	I3	4		AOHI60504T
				5		AOHI60505T
				7		AOHI60507T
	7.0	6.0	I4	4		AOHI70604T
				5		AOHI70605T
				7		AOHI70607T
	4.0	5.0	I1	4	Non-Hex	AOHI40504NT
				5		AOHI40505NT
				7		AOHI40507NT
	4.5	4.5	I2	4		AOHI45454NT
				5		AOHI45455NT
				7		AOHI45457NT
	6.0	5.0	I3	4		AOHI60504NT
				5		AOHI60505NT
				7		AOHI60507NT
	7.0	6.0	I4	4		AOHI70604NT
				5		AOHI70605NT
				7		AOHI70607NT



Type	MD (mm)	LB (mm)	Position No.	Height (mm)	Connection	Ref.C
Canine	5.0	6.5	C1	4	Regular	AOHC50654T
				5		AOHC50655T
				7		AOHC50657T
	5.0	6.5	C1	4	Non-Hex	AOHC50654NT
				5		AOHC50655NT
				7		AOHC50657NT




Type	MD (mm)	LB (mm)	Position No.	Height (mm)	Connection	Ref.C
Pre-Molar	4.5	6.0	P1	4	Regular	AOHM45604T
				5		AOHM45605T
				7		AOHM45607T
	5.0	7.0	P2	4		AOHM50704T
				5		AOHM50705T
				7		AOHM50707T
	4.5	6.0	P1	4	Non-Hex	AOHM45604NT
				5		AOHM45605NT
				7		AOHM45607NT
	5.0	7.0	P2	4		AOHM50704NT
				5		AOHM50705NT
				7		AOHM50707NT



Type	MD (mm)	LB (mm)	Position No.	Height (mm)	Connection	Ref. C
Molar	6.0	7.0	M1	4	Regular	AOHM60704T
				5		AOHM60705T
				7		AOHM60707T
	6.0	8.0	M2	4		AOHM60804T
				5		AOHM60805T
				7		AOHM60807T
	6.0	9.0	M3	4		AOHM60904T
				5		AOHM60905T
				7		AOHM60907T
	7.0	8.0	M4	4		AOHM70804T
				5		AOHM70805T
				7		AOHM70807T
	7.0	9.0	M5	4		AOHM70904T
				5		AOHM70905T
				7		AOHM70907T
	7.0	10.0	M6	4		AOHM70104T
				5		AOHM70105T
				7		AOHM70107T
	8.0	9.0	M7	4		AOHM80904T
				5		AOHM80905T
				7		AOHM80907T
	8.0	10.0	M8	4		AOHM80104T
				5		AOHM80105T
				7		AOHM80107T
	6.0	7.0	M1	4	Non-Hex	AOHM60704NT
				5		AOHM60705NT
				7		AOHM60707NT
	6.0	8.0	M2	4		AOHM60804NT
				5		AOHM60805NT
				7		AOHM60807NT
	6.0	9.0	M3	4		AOHM60904NT
				5		AOHM60905NT
				7		AOHM60907NT
	7.0	8.0	M4	4		AOHM70804NT
				5		AOHM70805NT
				7		AOHM70807NT
	7.0	9.0	M5	4		AOHM70904NT
				5		AOHM70905NT
				7		AOHM70907NT
	7.0	10.0	M6	4		AOHM70104NT
				5		AOHM70105NT
				7		AOHM70107NT
	8.0	9.0	M7	4		AOHM80904NT
				5		AOHM80905NT
				7		AOHM80907NT
	8.0	10.0	M8	4		AOHM80104NT
				5		AOHM80105NT
				7		AOHM80107NT
Special	4.5	6.0	S1	4	Regular	AOHS45604T
				5		AOHS45605T
				7		AOHS45607T
	5.0	6.5	S2	4		AOHS50654T
				5		AOHS50655T
				7		AOHS50657T
	5.0	7.0	S3	4		AOHS50704T
				5		AOHS50705T
				7		AOHS50707T
	6.0	7.0	S4	4		AOHS60704T
				5		AOHS60705T
				7		AOHS60707T
	6.0	8.0	S5	4		AOHS60804T
				5		AOHS60805T
				7		AOHS60807T
	6.0	9.0	S6	4		AOHS60904T
				5		AOHS60905T
				7		AOHS60907T
	7.0	8.0	S7	4		AOHS70804T
				5		AOHS70805T
				7		AOHS70807T
	7.0	9.0	S8	4		AOHS70904T
				5		AOHS70905T
				7		AOHS70907T
	7.0	10.0	S9	4		AOHS70104T
				5		AOHS70105T
				7		AOHS70107T
	8.0	9.0	S10	4		AOHS80904T
				5		AOHS80905T
				7		AOHS80907T
	8.0	10.0	S11	4		AOHS80104T
				5		AOHS80105T
				7		AOHS80107T
	4.5	6.0	S1	4	Non-Hex	AOHS45604NT
				5		AOHS45605NT
				7		AOHS45607NT
	5.0	6.5	S2	4		AOHS50654NT
				5		AOHS50655NT
				7		AOHS50657NT
	5.0	7.0	S3	4		AOHS50704NT
				5		AOHS50705NT
				7		AOHS50707NT
	6.0	7.0	S4	4		AOHS60704NT
				5		AOHS60705NT
				7		AOHS60707NT
	6.0	8.0	S5	4		AOHS60804NT
				5		AOHS60805NT
				7		AOHS60807NT
	6.0	9.0	S6	4		AOHS60904NT
				5		AOHS60905NT
				7		AOHS60907NT
	7.0	8.0	S7	4		AOHS70804NT
				5		AOHS70805NT
				7		AOHS70807NT
	7.0	9.0	S8	4		AOHS70904NT
				5		AOHS70905NT
				7		AOHS70907NT
	7.0	10.0	S9	4		AOHS70104NT
				5		AOHS70105NT
				7		AOHS70107NT
	8.0	9.0	S10	4		AOHS80904NT
				5		AOHS80905NT
				7		AOHS80907NT
	8.0	10.0	S11	4		AOHS80104NT
				5		AOHS80105NT
				7		AOHS80107NT

Abutment & Prosthetic Options

I. Fixture-level Prosthesis

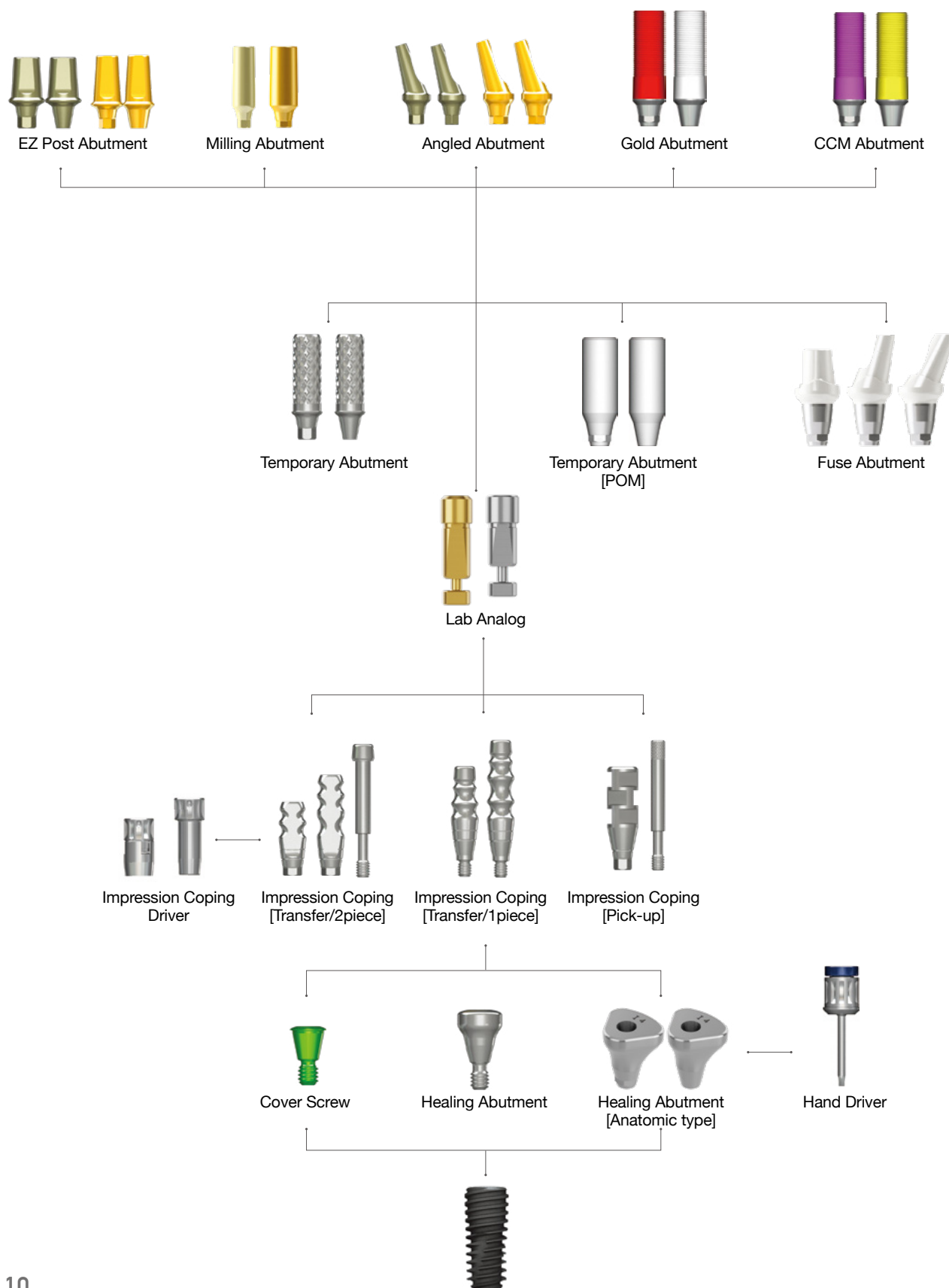


➡ Abutment Options & Components (Continued)

Impression Coping (2-piece, Transfer Type) (For closed-Tray Technique)

- Guide Pins : STTIG1611/ STTIG1615/ STTIG2011/
STTIG2015

- For use with Closed-tray technique Anti-rotation grooves are aligned with fixture hex direction
- Should be tightened with Impression Coping Driver
- Special impression coping screw which can be used with a 1.2mm hex driver is available on request.


Diameter	Type	Connection	Height (mm)	Ref.C
Ø4.0	2-piece	Mini	11	STTIM4011T
			15	STTIM4015T
Ø4.5			11	STTIM4511T
			15	STTIM4515T
Ø4.0		Regular	11	STTIS4011T
			15	STTIS4015T
Ø4.5			11	STTIS4511T
			15	STTIS4515T
Ø5.0			11	STTIS5011T
			15	STTIS5015T
Ø6.0			11	STTIS6011T
			15	STTIS6015T
Ø7.0			11	STTIS7011T
			15	STTIS7015T

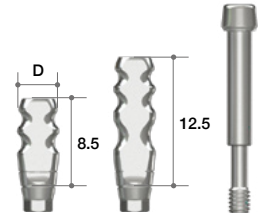

Impression Coping (1-piece, Transfer Type) (For closed-Tray Technique)

- Impression Coping design ensures easy and accurate transfer of fixture position.
- Impression Coping Driver and Hand Driver(1.2Hex) should be used to ensure Impression Coping is properly tightened.

Diameter	Type	Connection	Height (mm)	Ref.C
Ø4.0	1-piece	Mini	11	STTIM4011NH
Ø4.5			15	STTIM4015NH
			11	STTIM4511NH
			15	STTIM4515NH
		Regular	11	STTIS4011NH
15			STTIS4015NH	
Ø4.5			11	STTIS4511NH
			15	STTIS4515NH
Ø5.0			11	STTIS5011NH
			15	STTIS5015NH
Ø6.0			11	STTIS6011NH
			15	STTIS6015NH
Ø7.0			11	STTIS7011NH
			15	STTIS7015NH

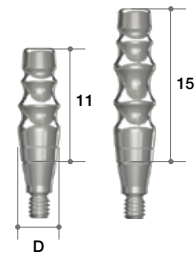


Impression Coping (2-piece, Pick-up Type) (For closed-Tray Technique)

- Guide Pins : STPIG1611(Mini / 11mm), STPIG1615(Mini/ 15mm),
STPIG2011(Regular/11mm),
STPIG2015(Regular/15mm) are included.

- Accurate impression-taking due to square structure
- Designed for easy & accurate impression-taking

Diameter	Type	Connection	Height(mm)	Ref.C		
Ø4.0	2-piece	Mini	Hex	11	STPIM4011T	
				15	STPIM4015T	
			Non-Hex	11	STPIM4011NT	
				15	STPIM4015NT	
Hex			11	STPIM4511T		
			15	STPIM4515T		
Non-Hex		11	STPIM4511NT			
		15	STPIM4515NT			
Ø4.0		Regular	Hex	11	STPIS4011T	
				15	STPIS4015T	
				Non-Hex	11	STPIS4011NT
					15	STPIS4015NT
Hex			11	STPIS4511T		
			15	STPIS4515T		
			Non-Hex	11	STPIS4511NT	
				15	STPIS4515NT	
Ø4.5	Hex		11	STPIS5011T		
			15	STPIS5015T		
			Non-Hex	11	STPIS5011NT	
				15	STPIS5015NT	
Ø5.0	Hex		11	STPIS6011T		
			15	STPIS6015T		
			Non-Hex	11	STPIS6011NT	
				15	STPIS6015NT	
Ø6.0	Hex		11	STPIS7011T		
			15	STPIS7015T		
			Non-Hex	11	STPIS7011NT	
				15	STPIS7015NT	
Ø7.0	Hex		11	STPIS7011T		
			15	STPIS7015T		
		Non-Hex	11	STPIS7011NT		
			15	STPIS7015NT		

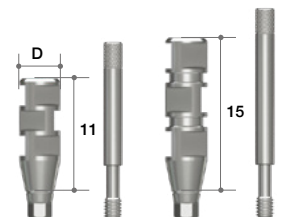


➡ Abutment Options & Components (Continued)

Bite Impression Coping

- Guide Pins: B1GP1624S/1634S/1644S/1646S/1666S/
2024S/2034S/2044S/2054S/2046S/2056S/2066S


- Bite Tray can be used to solve the opposing teeth and impressions at once
- By giving a groove for final tightening, it's possible to check the connection by X-ray
- Screw will not be tightened if the coping and the fixture are not properly connected
- Exclusive B.I.C Carrier is used to tight the Impression Coping conveniently
- Different color according to Cuff Size
- ※ Availability may vary by country

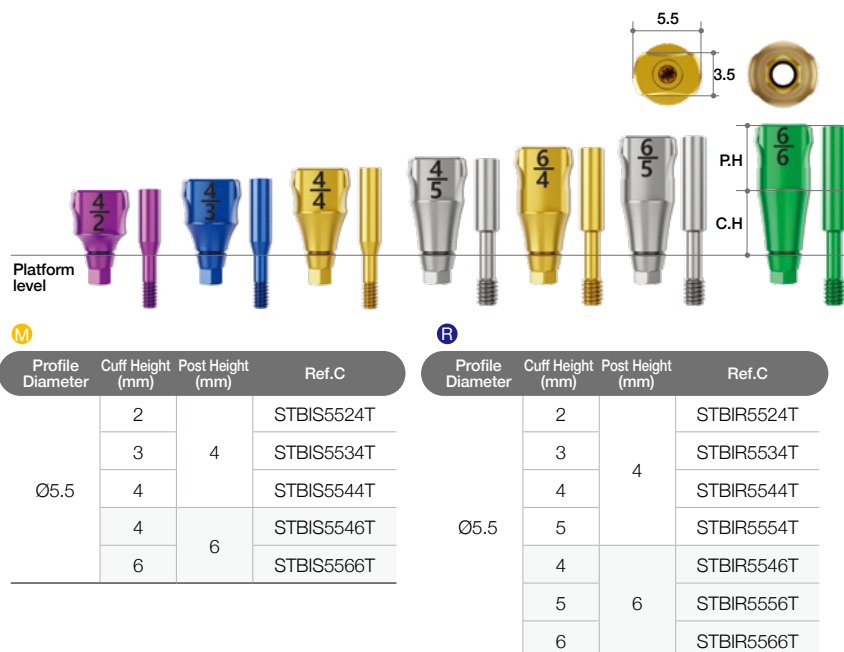

B.I.C Carrier

- Bite Impression Coping can be transferred and fastened at once
- ※ Availability may vary by country

Length (mm)	Ref.C
10	BICC10
14	BICC14

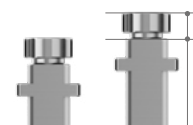


Analog

- Analog Screw(ALS18) included.

- For Chairside/ Labside
- Supporting Dental CAD
 - 3Shape
 - exocad
- 2 piece type



System		Profile Diameter	Height (mm)	Ref.C
ST	Mini	Ø4.0	10.5	STIALMT
	Regular	Ø4.0	9	STIALRT
Octa Level	Small	Ø3.8	9	OCTAALST
	Regular	Ø4.8		OCTAALRT
	Wide	Ø5.8		OCTAALWT
MUA Level (N-Type)		Ø4.8	9	MUAALT

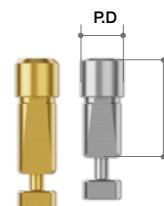


Temporary Abutment

- Abutment Screw is included
(Mini: STABSM / Regular: STABSS)

- To fabricate provisional restorations post area is designed to increase adhesive property of resin.
- Recommended torque: Mini 20Ncm/
Regular 25Ncm

Profile Diameter	Connection	Cuff Height(mm)	Ref.C
Ø4.0	Mini	Hex 1	STTA4010T
		3	STTA4030T
	Non-Hex	1	STTA4010NT
		3	STTA4030NT
Ø4.5	Regular	Hex 1	STTA4510T
		3	STTA4530T
	Non-Hex	1	STTA4510NT
		3	STTA4530NT

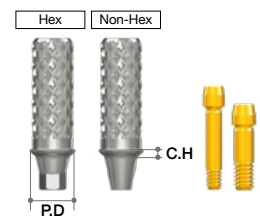


Temporary Abutment (POM)

- Abutment Screw(AS20) included

- For making chairside provisionals for the aesthetic zone.
- Especially useful for immediate placement after extraction.
- Available in both hex and non-hex.
- Recommend torque : 25Ncm

Profile Diameter	Cuff Height (mm)	Connection	Ref.C
Ø4.5	11.0	Hex	TA4511HPT
		Non-Hex	TA4511NPT

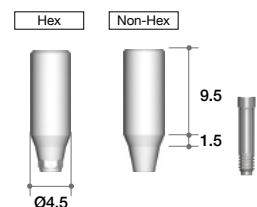


Fuse Abutment

- Abutment Screw(AS20)+Fuse cap included.

- Recommend torque : 25Ncm

Diameter	Cuff Height	Post Height	Type	Connection	Ref.C
Ø5.5	4	5.5	Straight	Regular	AOFAP5545P
	4	7	15°		AOFAP5415P
	4	7	25°		AOFAP5425P

NEW : 4mm cuff height available
→ Adequate for deeply placed implants
or thick gingival cases

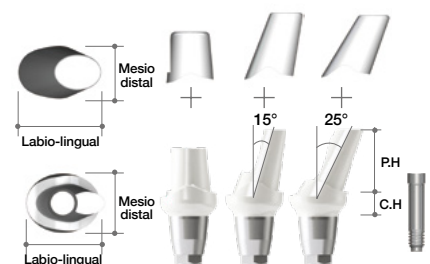

➔ Abutment Options & Components (Continued)


EZ Post Abutment

- Abutment Screw is included.

(Mini : STABSM / Regular : STABSS)

For cement-retained restoration:

- Post heights: 4.0, 5.5, 7.0mm
- Profile diameters: Ø4.5, Ø5.0, Ø6.0, Ø7.0
- Cuff heights: 1, 2, 3, 4, 5mm
- Aesthetic appearance realized via anodizing
- Non-hex abutments have no anti-rotation feature, and are contraindicated for single-unit restoration
- Recommended torque: Mini / 20Ncm, Regular / 30Ncm

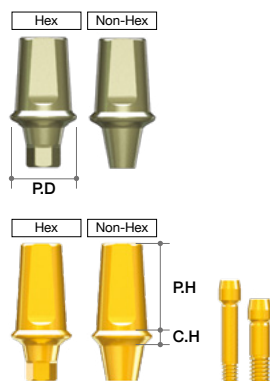


Profile Diameter	Connection	Cuff Height(mm)	Post Height(mm)	Ref.C	Profile Diameter	Connection	Cuff Height(mm)	Post Height(mm)	Ref.C
Ø4.5	Mini	Hex	1	STEA4611T	Ø5.0	Non-Hex	1	4.0	STEA5410NT
			2	STEA4621T			2	4.0	STEA5420NT
			3	STEA4631T			3	4.0	STEA5430NT
			4	STEA4641T			4	4.0	STEA5440NT
			5	STEA4651T			5	4.0	STEA5450NT
		Non-Hex	1	STEA4711T			1	5.5	STEA5610NT
			2	STEA4721T			2	5.5	STEA5620NT
			3	STEA4731T			3	5.5	STEA5630NT
			4	STEA4741T			4	5.5	STEA5640NT
			5	STEA4751T			5	5.5	STEA5650NT
	Regular	Hex	1	STEA4611NT		Non-Hex	1	7.0	STEA5710NT
			2	STEA4621NT			2	7.0	STEA5720NT
			3	STEA4631NT			3	7.0	STEA5730NT
			4	STEA4641NT			4	7.0	STEA5740NT
			5	STEA4651NT			5	7.0	STEA5750NT
Ø4.5	Regular	Hex	1	STEA4711NT	Ø6.0	Hex	1	4.0	STEA6410T
			2	STEA4721NT			2	4.0	STEA6420T
			3	STEA4731NT			3	4.0	STEA6430T
			4	STEA4741NT			4	4.0	STEA6440T
			5	STEA4751NT			5	4.0	STEA6450T
	Non-Hex	Hex	1	STEA4611T		Non-Hex	1	5.5	STEA6610T
			2	STEA4621T			2	5.5	STEA6620T
			3	STEA4631T			3	5.5	STEA6630T
			4	STEA4641T			4	5.5	STEA6640T
			5	STEA4651T			5	5.5	STEA6650T
	Regular	Hex	1	STEA4711T		Non-Hex	1	7.0	STEA6710T
			2	STEA4721T			2	7.0	STEA6720T
			3	STEA4731T			3	7.0	STEA6730T
			4	STEA4741T			4	7.0	STEA6740T
			5	STEA4751T			5	7.0	STEA6750T
	Non-Hex	Hex	1	STEA4611NT		Non-Hex	1	4.0	STEA6410NT
			2	STEA4621NT			2	4.0	STEA6420NT
			3	STEA4631NT			3	4.0	STEA6430NT
			4	STEA4641NT			4	4.0	STEA6440NT
			5	STEA4651NT			5	4.0	STEA6450NT
Ø5.0	Regular	Hex	1	STEA4711NT	Ø7.0	Hex	1	5.5	STEA6610NT
			2	STEA4721NT			2	5.5	STEA6620NT
			3	STEA4731NT			3	5.5	STEA6630NT
			4	STEA4741NT			4	5.5	STEA6640NT
			5	STEA4751NT			5	5.5	STEA6650NT
	Non-Hex	Hex	1	STEA4611T		Non-Hex	1	7.0	STEA6710NT
			2	STEA4621T			2	7.0	STEA6720NT
			3	STEA4631T			3	7.0	STEA6730NT
			4	STEA4641T			4	7.0	STEA6740NT
			5	STEA4651T			5	7.0	STEA6750NT
	Regular	Hex	1	STEA5410T		Non-Hex	1	5.5	STEA7610T
			2	STEA5420T			2	5.5	STEA7620T
			3	STEA5430T			3	5.5	STEA7630T
			4	STEA5440T			4	5.5	STEA7640T
			5	STEA5450T			5	5.5	STEA7650T
	Non-Hex	Hex	1	STEA5610T		Non-Hex	1	7.0	STEA7710T
			2	STEA5620T			2	7.0	STEA7720T
			3	STEA5630T			3	7.0	STEA7730T
			4	STEA5640T			4	7.0	STEA7740T
			5	STEA5650T			5	7.0	STEA7750T
Ø5.0	Regular	Hex	1	STEA5710T	Ø7.0	Non-Hex	1	5.5	STEA7610NT
			2	STEA5720T			2	5.5	STEA7620NT
			3	STEA5730T			3	5.5	STEA7630NT
			4	STEA5740T			4	5.5	STEA7640NT
			5	STEA5750T			5	5.5	STEA7650NT
	Non-Hex	Hex	1	STEA5410T		Non-Hex	1	7.0	STEA7710NT
			2	STEA5420T			2	7.0	STEA7720NT
			3	STEA5430T			3	7.0	STEA7730NT
			4	STEA5440T			4	7.0	STEA7740NT
			5	STEA5450T			5	7.0	STEA7750NT

Angled Abutment

- Abutment Screw is included
(Mini : STABSM / Regular : STABSS)

- Only one angle available: 17°
- Various diameters: Ø4.5, 5.0, 6.0
- Various cuff heights: 2mm, 4mm
- Covers 12 different angles (Hex, Hex-E)
- Aesthetic appearance realized via anodizing
- More convenient to use as milling area is widened by minimizing screw head height
- Recommended torque: Mini / 20Ncm, Regular / 30Ncm

Profile Diameter	Connection		Cuff Height(mm)	Post Height(mm)	Ref.C
Ø4.5	Mini	Hex-E	2	8	STAA4520MAT
			4		STAA4540MAT
		Hex	2		STAA4520MBT
			4		STAA4540MBT
		Non-Hex	2		STAA4520MNT
			4		STAA4540MNT
Ø4.5	Regular	Hex-E	2	8	STAA4520AT
			4		STAA4540AT
		Hex	2		STAA4520BT
			4		STAA4540BT
		Non-Hex	2		STAA4520NT
			4		STAA4540NT
Ø5.0	Regular	Hex-E	2		STAA5020AT
			4		STAA5040AT
		Hex	2		STAA5020BT
			4		STAA5040BT
		Non-Hex	2		STAA5020NT
			4		STAA5040NT
Ø6.0	Regular	Hex-E	2		STAA6020AT
			4		STAA6040AT
		Hex	2		STAA6020BT
			4		STAA6040BT
		Non-Hex	2		STAA6020NT
			4		STAA6040NT

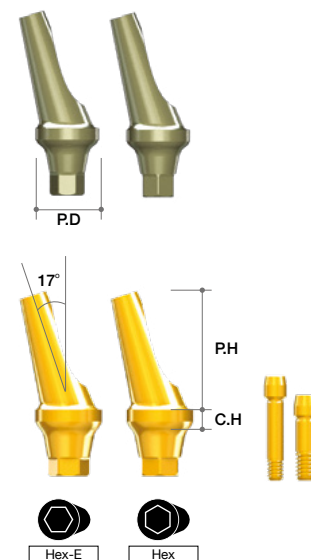


Milling Abutment

- Abutment Screw is included
(Mini : STABSM / Regular : STABSS)

- Used for customized milling
- Profile diameters: Ø4.5/ Ø5.0/ Ø6.0/ Ø7.0
- Cuff heights: 1.5mm, 3.0mm
- Aesthetic appearance realized via anodizing
- Recommended torque: Mini / 20Ncm, Regular / 30Ncm

Profile Diameter	Connection	Cuff Height(mm)	Post Height(mm)	Ref.C
Ø4.5	Mini	1.5	10.5	STMA4015T
		3.0	9.0	STMA4030T
Ø4.5	Regular	1.5	10.5	STMAS4015T
		3.0	9.0	STMAS4030T
Ø5.0	Regular	1.5	10.5	STMAS5015T
		3.0	9.0	STMAS5030T
Ø5.0	Regular	1.5	8.0	STMASG5015T
		3.0	7.0	STMASG5030T
Ø6.0	Regular	1.5	10.5	STMASG6015T
		3.0	9.0	STMASG6030T
Ø7.0	Regular	1.5	10.5	STMASG7015T
		3.0	9.0	STMASG7030T

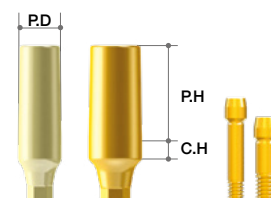


➔ Abutment Options & Components

Milling Abutment Type II (GALLI Abutment)

- AnyOne Internal : Abutment Screw (AS20) included.

- Long post enables easier customization from milling.
- Recommend torque : 35Ncm

Profile Diameter	Connection	Cuff Height (mm)	Post Height (mm)	Ref.C
Ø4.0	Regular	1.0	19	AOBOT4019HT
				AOBOT4019NT

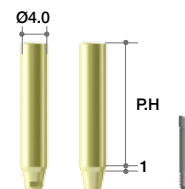


Gold Abutment

- Abutment Screw(AS20) included

- For fabrication of customized abutment for either screw or cement retained restorations.
- Available in both hex (red) and non-hex (white)
- Melting point of gold alloy : 1063°C
- Threaded sleeves allow for better retention of resin or wax.
- Recommend torque : 30Ncm

Profile Diameter	Connection	Cuff Height (mm)	Post Height (mm)	Ref.C
Ø4.5	Regular	1.0	11.0	GA4515HT
				GA4515NT

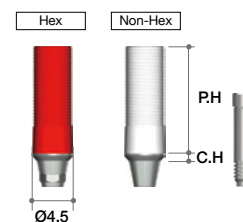


CCM Abutment

- Abutment Screw(AS20) included


- Useful to make a customized abutment in difficult situations.
- Can be casted with non-precious alloys(Ni-Cr, Cr-Co alloys).
- Non-precious melting temperature : Depend on Manufacturer
- Threaded sleeves for convenient Resin / Wax-up.
- Melting temperature of CCM : 1300~1400°C
- Recommend torque : 35Ncm

Profile Diameter	Connection	Cuff Height (mm)	Post Height (mm)	Ref.C
Ø4.5	Regular	1.0	11.0	CA4515HT
				CA4515NT

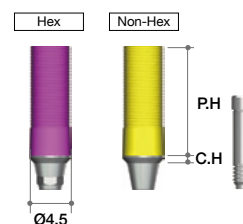

I. Fixture Level Prosthesis

1. Fixture Level Prosthesis_Digital

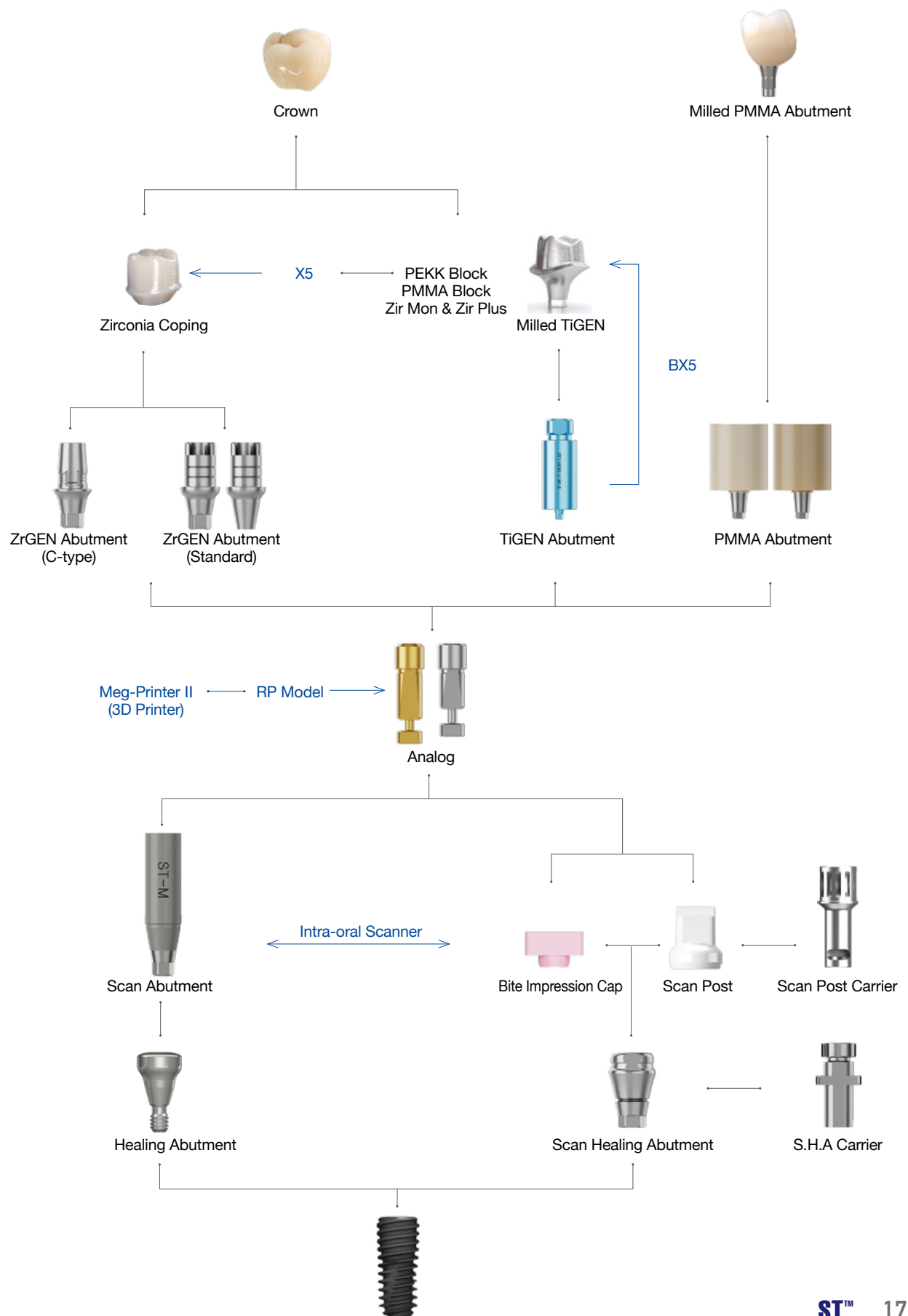


» Digital Material

I. ZrGEN®

ZrGEN® is the brand name of MegaGen's Titanium Base. ZrGEN provides an aesthetic outcome and simplified dental implant prosthesis. A ZrGEN® crown and monolithic crown connected to a ZrGEN® Abutment provide strong and precise connection with the implant fixture.

Variety of ZrGEN®

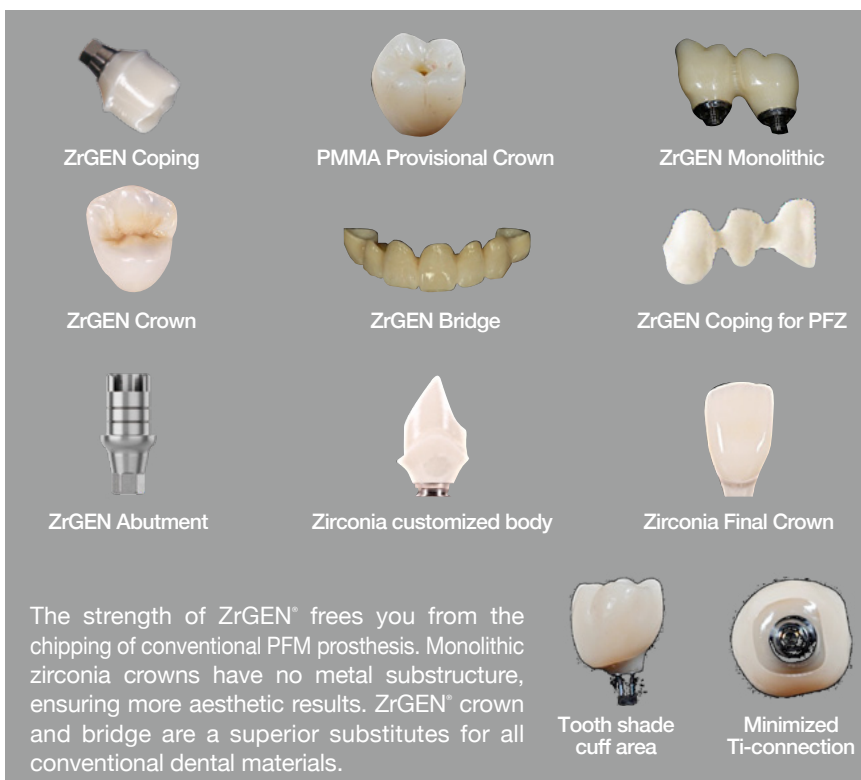

ZrGEN® Sub Structure

ZrGEN®

II. TiGEN®

TiGEN® is the brand name of MegaGen's Pre-milled Abutment. It promises outstanding durability and simplified dental implant prosthesis. Ready-made connection part provides a strong and precise connection with the implant fixture.



Clinical Application




➡ Abutment & Prosthetic Option (Continued)

Scan Abutment

- Abutment Screw included.
 - ✓ ST Mini (SAOSGSAS3110)
 - ✓ ST Regular (SAOSGSAS3210)
 - ✓ Octa Level (SAIRCS200)
 - ✓ MUA Level (SAMUAS)
- For Chairside/ Labside
- Included spare Abutment Screw
- Supporting Dental CAD
 - 3Shape
 - exocad
 - Dental Wings
- Recommend torque : By Hand (5~8Ncm)

System		Profile Diameter	Height (mm)	Ref.C
ST	Mini	Ø4.0	9	OSGSSC3110T
			13	OSGSSC3111T
	Regular		9	OSGSSC3210T
			13	OSGSSC3211T
Octa Level		Ø4.0	11	AOCESC4011T
MUA Level (N Type)		Ø4.0	13	AMUASR4013T

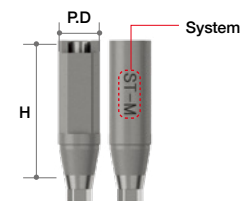

Scan Healing Abutment

- Abutment Screw included.
 - ✓ ST Mini (STMHS1604/ STMHS1605/ ST-MHS1607/ STMHS1609)
 - ✓ ST Regular (STRHS2004/ STRHS2005/ STRHS2007/ STRHS2009)
- Scannable Healing Abutment.
- For accurate scanning, Scan Healing Abutment must be exposed at least 2.0mm from surgical site.
- Profile Diameter can be checked by the number of Groove.
 - Profile Diameter : Ø4 → Groove : 0ea
 - Profile Diameter : Ø5 → Groove : 1ea
 - Profile Diameter : Ø6 → Groove : 2ea
 - Profile Diameter : Ø7 → Groove : 3ea
- Height can be checked by the number of Notch.
 - Height : 4mm → Notch : 0ea
 - Height : 5mm → Notch : 1ea
 - Height : 7mm → Notch : 2ea
 - Height : 9mm → Notch : 3ea
- Recommend torque : By Hand (5~8Ncm)
- Height 9mm & ST system - FDA : Approved in 2023



System	Profile Diameter	Height (mm)	Ref.C
ST	Mini	Ø4.0	4 STMHS4004T
		5	STMHS4005T
		7	STMHS4007T
		9	STMHS4009T
	Regular	Ø5.0	4 STRSH5004T
		5	STRSH5005T
		7	STRSH5007T
		9	STRSH5009T
		Ø6.0	4 STRSH6004T
		5	STRSH6005T
		7	STRSH6007T
		9	STRSH6009T
	Regular	Ø7.0	4 STRSH7004T
		5	STRSH7005T
		7	STRSH7007T
		9	STRSH7009T

S.H.A Carrier

(Coming Soon)

- It is used by fastening it to the head of the scan healing abutment.
- It is used when transporting the S.H.A Carrier to the fixture after fastening it to the Scan Healing Abutment.

Diameter	Length (mm)	Ref.C
Ø4.0	10	SHC4010
	14	SHC4014
Ø5.0	10	SHC5010
	14	SHC5014
Ø6.0	10	SHC6010
	14	SHC6014
Ø7.0	10	SHC7010
	14	SHC7014




➡ Abutment & Prosthetic Option

Scan Post

- Scan Healing Abutment should be exposed 2.0mm from the surgical site for accurate scanning. Scanning would be much easier if you connect Scan Post when scanning seems difficult due to less exposure of Scan Healing Abutment or other conditions.
- Select Scan Post based on the diameter of Scan Healing Abutment
- Scan Post is a disposable product.
- 1 set consists of 10 Scan Posts.


Profile Diameter	Height (mm)	Ref.C
Ø4.0	6.5	SP4007.MTN
Ø5.0		SP5007.MTN
Ø6.0		SP6007.MTN
Ø7.0		SP7007.MTN

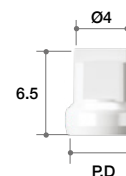

Scan Post Carrier

- It is used by fastening it to the head of the scan post.
- It is used when transporting the scan post to the fixture after fastening it to the scan healing abutment.



Profile Diameter	Height (mm)	Ref.C
Ø5	19	SPC16

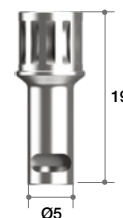

Bite Impression Cap

- If the Scan Healing Abutment is not exposed to more than 2.0mm at the surgical site or is not easy to scan due to environment such as interference of the surrounding values, the Bite Impression Cap is fastened to obtain the Scan, Impression, or Bite.
- Bite Impression Cap is selected according to the Profile Diameter of the Scan Healing Abutment.
- The Bite Impression Cap is for one-time use and includes 5 in 1 set.
- ※ Availability may vary by country

Profile Diameter	Height (mm)	Ref.C
Ø4.0	2	BIC4002P
	4	BIC4004P
	6	BIC4006P
Ø5.0	2	BIC5002P
	4	BIC5004P
	6	BIC5006P
Ø6.0	2	BIC6002P
	4	BIC6004P
	6	BIC6006P
Ø7.0	2	BIC7002P
	4	BIC7004P
	6	BIC7006P

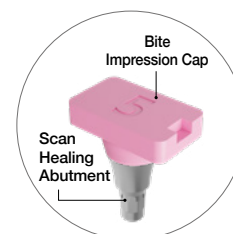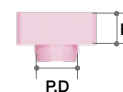


➡ Analog Option

Analog

- Analog Screw(ALS18) included.
- For Chairside/ Labside
- Supporting Dental CAD
 - 3Shape
 - exocad
- 2 piece type



System	Profile Diameter	Height (mm)	Ref.C
ST	Mini	Ø4.0 10.5	STIALMT
	Regular	Ø4.0 9	STIALRT

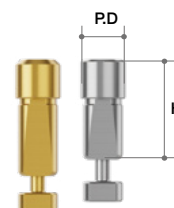


➡ PMMA Abutment Option

PMMA Abutment

(Coming soon)

- Abutment Screw included.
 - ✓ ST Mini (OSGSAS3110)
 - ✓ ST Regular (OSGSAS3210)
- Regular : Compatible with AnyOne Internal system
- Pre-milled Abutment
 - Pre-milled part : Implant Connection + Cuff (0.6/ 1.5/ 3.0mm)
- 1 Set consists of 10 Abutments
 - included spare Abutment Screw
- Supporting Dental CAD
 - 3 Shape
 - exocad
- Supporting Milling Machine
 - MegaGen Implant : BX5
 - ARUM DENTISTRY
- Recommend torque : 20Ncm
- FDA : Approved in 2024
- CE : Approved in 2024

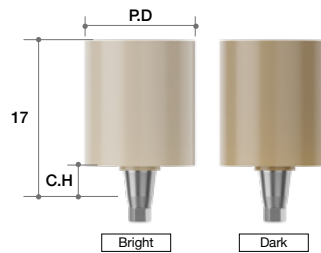


[PMMA Abutment Line-Up]

	Bright		
	D8	D12	D16
C0.6			
C1.5			
C3.0			
	Dark		
	D8	D12	D16
C0.6			
C1.5			
C3.0			

System	Color	Profile Diameter	Cuff Height (mm)	Length (mm)	Type	Ref.C
ST	Bright	Ø8	0.6	17	Hex	STPAM0608B.MTN
			1.5			STPAM1508B.MTN
			3.0			STPAM3008B.MTN
		Ø12	0.6			STPAM0612B.MTN
			1.5			STPAM1512B.MTN
			3.0			STPAM3012B.MTN
		Ø16	0.6			STPAM0616B.MTN
			1.5			STPAM1516B.MTN
			3.0			STPAM3016B.MTN
	Dark	Ø8	0.6			STPAM0608D.MTN
			1.5			STPAM1508D.MTN
			3.0			STPAM3008D.MTN
		Ø12	0.6			STPAM0612D.MTN
			1.5			STPAM1512D.MTN
			3.0			STPAM3012D.MTN
		Ø16	0.6			STPAM0616D.MTN
			1.5			STPAM1516D.MTN
			3.0			STPAM3016D.MTN
		Ø8	0.6	17	Hex	AOPA0608B.MTN
			1.5			AOPA1508B.MTN
			3.0			AOPA3008B.MTN
		Ø12	0.6			AOPA0612B.MTN
			1.5			AOPA1512B.MTN
			3.0			AOPA3012B.MTN
		Ø16	0.6			AOPA0616B.MTN
			1.5			AOPA1516B.MTN
			3.0			AOPA3016B.MTN
	Regular	Ø8	0.6			AOPA0608D.MTN
			1.5			AOPA1508D.MTN
			3.0			AOPA3008D.MTN
		Ø12	0.6			AOPA0612D.MTN
			1.5			AOPA1512D.MTN
			3.0			AOPA3012D.MTN
		Ø16	0.6			AOPA0616D.MTN
			1.5			AOPA1516D.MTN
			3.0			AOPA3016D.MTN

➡ ZrGEN Abutment Option (Continued)

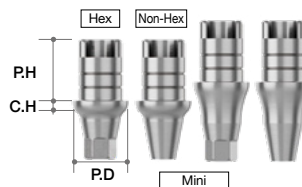

ZrGEN Abutment

- Abutment Screw included.

✓ ST Mini (OSGSAS3110)

✓ ST Regular (OSGSAS3210)


✓ Octa Level (IRCS200)

✓ MUA Level (MUAS)


- Titanium Base
- 1 set consists of 10 Abutments.
 - included spare Abutment Screw
 - Mini ZrGEN has special abutment screw (Available only in ZrGEN Abutment)
- Supporting DentalCAD
 - 3 Shape
 - exocad
- - Dental Wing
- Post Height can be checked by the number of Groove.
 - Post Height : 4.5mm → Groove : 2ea
 - Post Height : 5mm → Groove : 3ea
 - Post Height : 6mm → Groove : 4ea
 - Post Height : 8mm → Groove : 6ea
- Recommend torque
 - 30Ncm : ST (Regular)
 - 20Ncm : ST (Mini)
 - 35Ncm : Octa Level
 - 15Ncm : MUA Level

Standard

System	Profile Diameter	Cuff Height (mm)	Post Height (mm)	Type	Ref.C
ST	Mini	Ø4.0	4.5	Hex	OSGSPA3111.MTN
			6		OSGSPA3112.MTN
			8		OSGSPA3113.MTN
			4.5		OSGSPA3121.MTN
			6		OSGSPA3122.MTN
			8		OSGSPA3123.MTN
			4.5		OSGSPA3131.MTN
			6		OSGSPA3132.MTN
			8		OSGSPA3133.MTN
			4.5		OSGSPA3141.MTN
			6		OSGSPA3142.MTN
			8		OSGSPA3143.MTN
			4.5	Non-Hex	OSGSPA3111N.MTN
			6		OSGSPA3112N.MTN
			8		OSGSPA3113N.MTN
			4.5		OSGSPA3121N.MTN
			6		OSGSPA3122N.MTN
			8		OSGSPA3123N.MTN
			4.5		OSGSPA3131N.MTN
			6		OSGSPA3132N.MTN
			8		OSGSPA3133N.MTN
			4.5		OSGSPA3141N.MTN
			6		OSGSPA3142N.MTN
			8		OSGSPA3143N.MTN

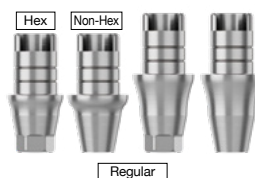

System	Profile Diameter	Cuff Height (mm)	Post Height (mm)	Type	Ref.C
Regular	Ø4.0	0.6	4.5	Hex	OSGSPA3211.MTN
			6		OSGSPA3212.MTN
			8		OSGSPA4018.MTN
		1.5	4.5		OSGSPA4025.MTN
			6		OSGSPA4026.MTN
			8		OSGSPA4028.MTN
		3.0	4.5		OSGSPA4035.MTN
			6		OSGSPA4036.MTN
			8		OSGSPA4038.MTN
		4.0	4.5		OSGSPA4045.MTN
			6		OSGSPA4046.MTN
			8		OSGSPA4048.MTN
		0.6	4.5	Non-Hex	OSGSPA3211N.MTN
			6		OSGSPA3212N.MTN
			8		OSGSPA4018N.MTN
		1.5	4.5		OSGSPA4025N.MTN
			6		OSGSPA4026N.MTN
			8		OSGSPA4028N.MTN
		3.0	4.5		OSGSPA4035N.MTN
			6		OSGSPA4036N.MTN
			8		OSGSPA4038N.MTN
		4.0	4.5		OSGSPA4045N.MTN
			6		OSGSPA4046N.MTN
			8		OSGSPA4048N.MTN
ST	Regular	0.6	4.5	Hex	OSGSPA4515.MTN
			6		OSGSPA4516.MTN
			8		OSGSPA4518.MTN
		1.5	4.5		OSGSPA3221.MTN
			6		OSGSPA3222.MTN
			8		OSGSPA4528.MTN
		3.0	4.5		OSGSPA4535.MTN
			6		OSGSPA4536.MTN
			8		OSGSPA4538.MTN
		4.0	4.5		OSGSPA4545.MTN
			6		OSGSPA4546.MTN
			8		OSGSPA4548.MTN
		0.6	4.5	Non-Hex	OSGSPA4515N.MTN
			6		OSGSPA4516N.MTN
			8		OSGSPA4518N.MTN
		1.5	4.5		OSGSPA3221N.MTN
			6		OSGSPA3222N.MTN
			8		OSGSPA4528N.MTN
		3.0	4.5		OSGSPA4535N.MTN
			6		OSGSPA4536N.MTN
			8		OSGSPA4538N.MTN
		4.0	4.5		OSGSPA4545N.MTN
			6		OSGSPA4546N.MTN
			8		OSGSPA4548N.MTN

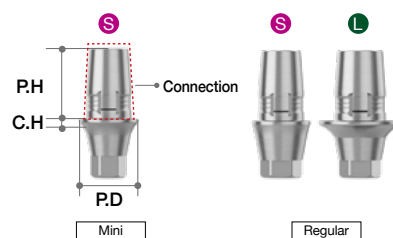

- Abutment Screw included.
- ✓ ST Mini (OSGSAS3110)
- ✓ ST Regular (OSGSAS3210)

- Titanium base for CEREC users.
- In in Lab CAD Software, compatible with Xive Library.
- 1 set consists of 10 Abutments.
 - included spare Abutment Screw.
- Recommend torque
 - 30Ncm : ST (Regular)
 - 20Ncm : ST (Mini)



C-type

System	Profile Diameter	Cuff Height (mm)	Post Height	Connection	Ref.C
ST	Mini	Ø3.9	0.5	S	STCSS3405.MTN
			1		STCSS3410.MTN
			2		STCSS3420.MTN
		Ø4.3	0.5		STCSS3805.MTN
			1		STCSS3810.MTN
			2		STCSS3820.MTN
	Regular	Ø3.9	0.5	S	STCSR3405.MTN
			1		STCSR3410.MTN
			2		STCSR3420.MTN
		Ø4.3	0.5		STCSR3805.MTN
			1		STCSR3810.MTN
			2		STCSR3820.MTN
		Ø5.5	0.5	L	STCLR4505.MTN
			1		STCLR4510.MTN
			2		STCLR4520.MTN


➡ TiGEN Abutment Option (Continued)

TiGEN Abutment

- Abutment Screw included.
 - ✓ ST Mini (OSGSAS3110)
 - ✓ ST Regular (OSGSAS3210)
- Pre-milled Abutment
- 1 set consists of 10 Abutments.
 - included spare Abutment Screw
- Supporting DentalCAD
 - 3 Shape
 - exocad
 - Dental Wings
- Supporting Milling Machine
 - MegaGen Implant : BX5
 - ARUM DENTISTRY
- Recommend torque
 - 30Ncm : ST (Regular)
 - 20Ncm : ST (Mini)
- FDA : Approved in 2023

MegaGen type


System		Color	Profile Diameter	Height (mm)	Type	Ref.C		
ST	Mini	Sky blue	Ø10	20	Hex	OSTG3112.MTN		
					Non-Hex	OSTG3112N.MTN		
					Hex	OSTG3111.MTN		
					Non-Hex	OSTG3111N.MTN		
	Regular				Ø10	Hex	OSTG3212.MTN	
						Non-Hex	OSTG3212N.MTN	
						Ø12	Hex	OSTG3211.MTN
							Non-Hex	OSTG3211N.MTN

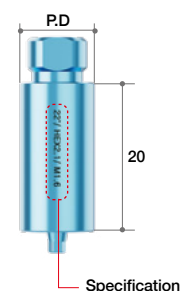

- Abutment Screw included.
 - ✓ ST Mini (OSGSAS3110)
 - ✓ ST Regular (OSGSAS3210)
- Pre-milled Abutment
- 1 set consists of 10 Abutments.
 - included spare Abutment Screw
- Supporting DentalCAD
 - 3 Shape
 - exocad
- Recommend torque
 - 30Ncm : ST (Regular)
 - 20Ncm : ST (Mini)
- FDA : Approved in 2023

NT type


System		Color	Profile Diameter	Height (mm)	Type	Ref.C
ST	Mini	Sky blue	Ø10	16	Hex	OSTGN3112.MTN
					Non-Hex	OSTGN3112N.MTN
			Ø12		Hex	OSTGN3111.MTN
					Non-Hex	OSTGN3111N.MTN
	Regular		Ø10		Hex	OSTGN3212.MTN
					Non-Hex	OSTGN3212N.MTN
			Ø12		Hex	OSTGN3211.MTN
					Non-Hex	OSTGN3211N.MTN

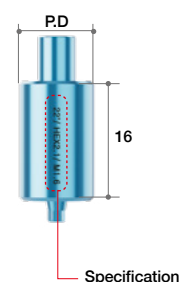

- Abutment Screw included.
 - ✓ ST Mini (OSGSAS3110)
 - ✓ ST Regular (OSGSAS3210)

Medentika type

System	Color	Profile Diameter	Height (mm)	Type	Ref.C
ST	Mini	Sky blue	Ø12	Hex	OSTGM3111.MTN
				Non-Hex	OSTGM3111N.MTN
				Hex	OSTGM3211.MTN
				Non-Hex	OSTGM3211N.MTN

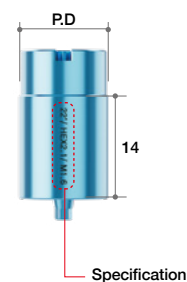


➡ TiGEN Abutment Option

- Abutment Screw included.
- ✓ ST Mini (OSGSAS3110)
- ✓ ST Regular (OSGSAS3210)

- Pre-milled Abutment
 - Pre-milled part : Implant Connection + Cuff (0.6/ 1.5/ 3.0mm)
- 1 set consists of 10 Abutments
 - included spare Abutment Screw
- Used by fastening it to a reverse jig
- Supporting Dental CAD
 - 3 Shape
 - exocad
- Supporting Milling Machine
 - MegaGen Implant : BX5
 - ARUM DENTISTRY
- Recommend torque
 - 30Ncm : ST (Regular)
 - 20Ncm : ST (Mini)
- FDA : Approved in 2023
- CE : Approved in 2024

CUFF type

System	Color	Profile Diameter	Cuff Height (mm)	Height (mm)	Type	Ref.C
ST	Mini	Ø8	0.6	16	Hex	OSTGMR0608.MTN
			1.5			OSTGMR1508.MTN
			3.0			OSTGMR3008.MTN
		Ø10	0.6			OSTGMR0610.MTN
			1.5			OSTGMR1510.MTN
			3.0			OSTGMR3010.MTN
		Ø12	0.6			OSTGMR0612.MTN
			1.5			OSTGMR1512.MTN
			3.0			OSTGMR3012.MTN
	Regular	Ø8	0.6			OSTGRR0608.MTN
			1.5			OSTGRR1508.MTN
			3.0			OSTGRR3008.MTN
		Ø10	0.6			OSTGRR0610.MTN
			1.5			OSTGRR1510.MTN
			3.0			OSTGRR3010.MTN
		Ø12	0.6			OSTGRR0612.MTN
			1.5			OSTGRR1512.MTN
			3.0			OSTGRR3012.MTN

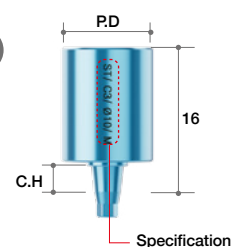

[TiGEN Abutment CUFF type Line-Up]





	D8	D10	D12
C0.6			
C1.5			
C3.0			

➤ CUFF types of TiGEN Abutment have same form of cuff shape as the Scan Healing Abutment thus custom abutment with perfectly fit to emergence profile can be fabricated

· Various cuff sizes for every gingiva height

➤ Pre-milled cuff reduces milling time + precision is increased with reverse jig milling

- 60% reduction in milling time when compared with conventional products!
- NO post milling, allowing reverse jig milling to occlusal surface within 8 minutes!

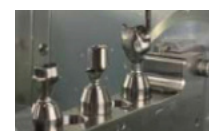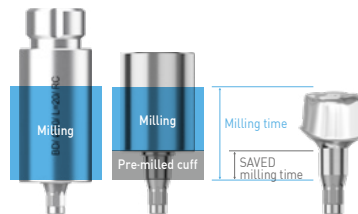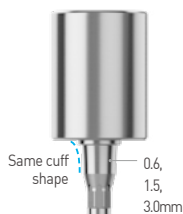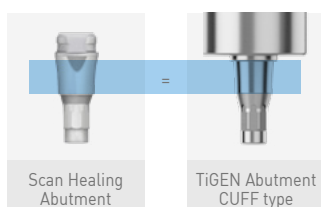

➡ Reverse Jig Connector Option

Reverse Jig Connector

- Milling screws exclusively for Reverse Jig Connector are included


- ✓ ST Mini (STRJMSM)
- ✓ ST Regular (STRJMSR)

- Do not use Non-Engage(Hex/ Octa)
- System can be checked by the number of **Groove**
 - ST Mini → Groove : narrow 3ea
 - ST Regular → Groove : wide 3ea
- Available milling machines
 - MegaGen Implant : BX5
 - ARUM DENTISTRY
- Recommended Torque
 - 35Ncm
 - Dedicated Driver (DP-RV-TORQ-DRV) (option)
- When Connected counterpart to Reverse Jig use Allen Wrench
 - Allen Key Size : 2.5mm
 - Dedicated Wrench (DP-HEX-TWLENCH) (option)


System	Type	Ref.C
ST	Mini	STTGRJM00P
	Regular	STTGRJR00P



*An example of grooves for different systems

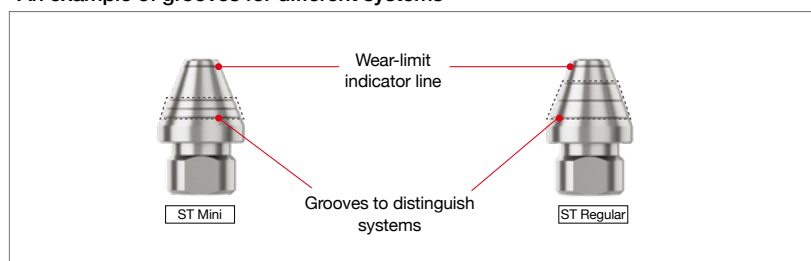


*Reverse Jig Connector assembly Option

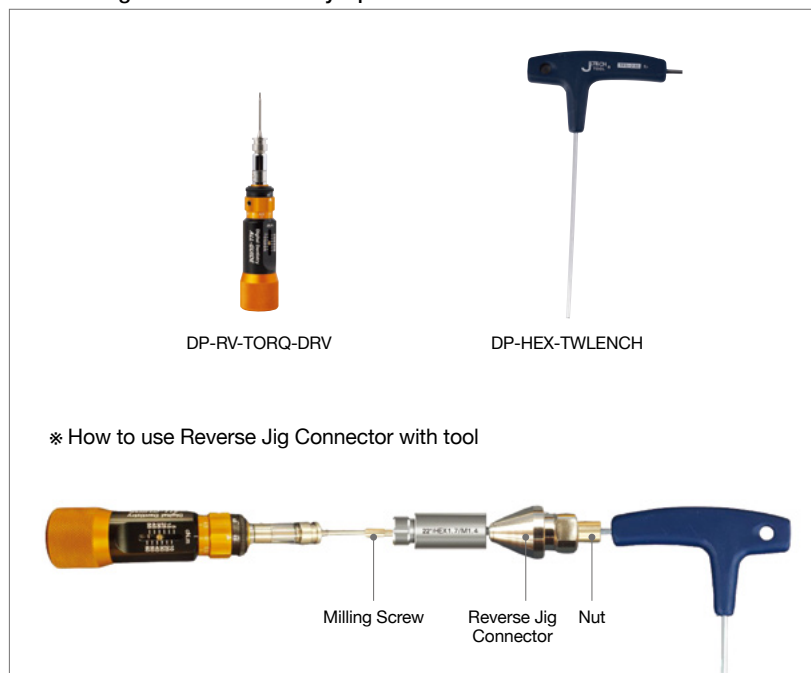

II. Abutment-level Prosthesis

1. Solid Abutment & Components

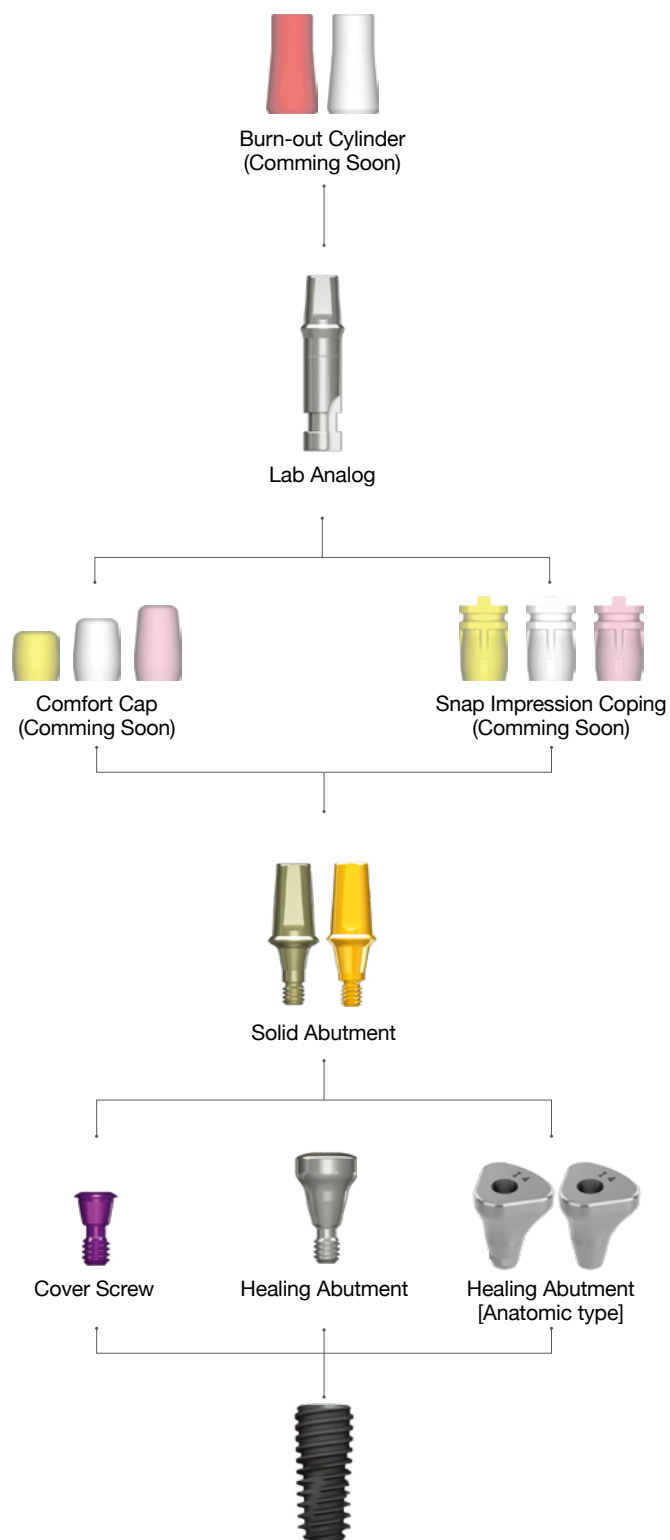

➔ Abutment Options & Components

Solid Abutment

For cement-retained restoration:

- Impression should be taken after connecting abutment to fixture
- Connect using 1.2mm hand driver
- Recommended torque: Mini / 20Ncm, Regular / 30Ncm
- Profile diameters: Ø4.0, Ø4.5, Ø5.0, Ø6.0, Ø7.0
- Cuff heights: 1.0, 2.0, 3.0, 4.0, 5.0mm
- Post heights: 4.0, 5.5, 7.0mm



Profile Diameter	Connection	Cuff Height(mm)	Post Height(mm)	Ref.C
Ø4.0	Mini	1	4.0	STSA4410
		2		STSA4420
		3		STSA4430
		4		STSA4440
		5		STSA4450
		1	5.5	STSA4610
		2		STSA4620
		3		STSA4630
		4		STSA4640
		5		STSA4650
		1	7.0	STSA4710
		2		STSA4720
		3		STSA4730
		4		STSA4740
		5		STSA4750
Ø4.5	Mini	1	4.0	STSA4411
		2		STSA4421
		3		STSA4431
		4		STSA4441
		5		STSA4451
		1	5.5	STSA4611
		2		STSA4621
		3		STSA4631
		4		STSA4641
		5		STSA4651
		1	7.0	STSA4711
		2		STSA4721
		3		STSA4731
		4		STSA4741
		5		STSA4751
Profile Diameter	Connection	Cuff Height(mm)	Post Height(mm)	Ref.C
Ø4.0	Regular	1	4.0	STSAS4410
		2		STSAS4420
		3		STSAS4430
		4		STSAS4440
		5		STSAS4450
		1	5.5	STSAS4610
		2		STSAS4620
		3		STSAS4630
		4		STSAS4640
		5		STSAS4650
		1	7.0	STSAS4710
		2		STSAS4720
		3		STSAS4730
		4		STSAS4740
		5		STSAS4750
Ø4.5	Regular	1	4.0	STSAS4411
		2		STSAS4421
		3		STSAS4431
		4		STSAS4441
		5		STSAS4451
		1	5.5	STSAS4611
		2		STSAS4621
		3		STSAS4631
		4		STSAS4641
		5		STSAS4651
		1	7.0	STSAS4711
		2		STSAS4721
		3		STSAS4731
		4		STSAS4741
		5		STSAS4751
Profile Diameter	Connection	Cuff Height(mm)	Post Height(mm)	Ref.C
Ø5.0	Regular	1	4.0	STSA5410
		2		STSA5420
		3		STSA5430
		4		STSA5440
		5		STSA5450
		1	5.5	STSA5610
		2		STSA5620
		3		STSA5630
		4		STSA5640
		5		STSA5650
		1	7.0	STSA5710
		2		STSA5720
		3		STSA5730
		4		STSA5740
		5		STSA5750
	Regular	1	4.0	STSA6410
		2		STSA6420
		3		STSA6430
		4		STSA6440
		5		STSA6450
		1	5.5	STSA6610
		2		STSA6620
		3		STSA6630
		4		STSA6640
		5		STSA6650
Ø6.0	Regular	1	5.5	STSA6610
		2		STSA6620
		3		STSA6630
		4		STSA6640
		5		STSA6650
		1	7.0	STSA6710
		2		STSA6720
		3		STSA6730
		4		STSA6740
		5		STSA6750
Ø7.0	Regular	1	5.5	STSA7610
		2		STSA7620
		3		STSA7630
		4		STSA7640
		5		STSA7650


Lab Analog

Diameter	Height(mm)	Ref.C
Ø4.0	4.0	STRLA440
	5.5	STRLA460
	7.0	STRLA470
Ø4.5	4.0	STRLA441
	5.5	STRLA461
	7.0	STRLA471
Ø5.0	4.0	STRLA540
	5.5	STRLA560
	7.0	STRLA570
Ø6.0	4.0	STRLA640
	5.5	STRLA660
	7.0	STRLA670
Ø7.0	5.5	STRLA760



II. Abutment Level Prosthesis

2. Octa Abutment & Components

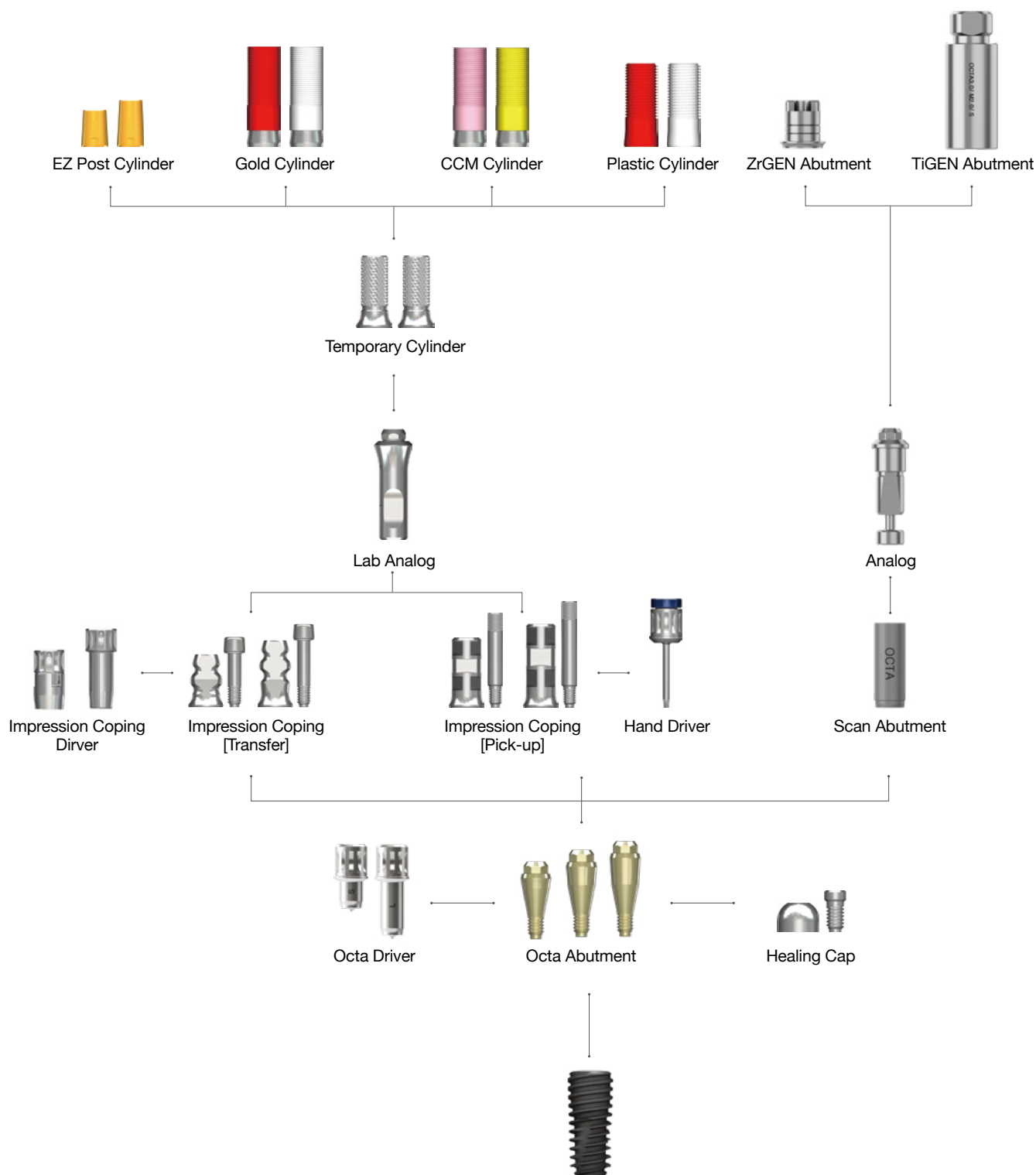

➔ Components for Octa Abutment (Continued)

Octa Abutment

- Used to make multiple screw-retained prosthetics.
- Recommend torque : 35Ncm

Profile Diameter	Connection	Cuff Height(mm)	Ref.C
Ø3.8	Regular	1.0	OA4010
		1.5	OA4015
		2.5	OA4025
		3.5	OA4035
		4.5	OA4045
5.5		OA4055	
Ø4.8		1.0	OA5010
		1.5	OA5015
		2.5	OA5025
		3.5	OA5035
		4.5	OA5045
Ø5.8		5.5	OA5055
		1.0	OA6010
		1.5	OA6015
		2.5	OA6025
		3.5	OA6035
		4.5	OA6045
		5.5	OA6055

P.D

C.H

P.D

P.D

Healing Cap

- Cylinder Screw (IRCS200) included
- Protects Octa Abutment and minimizes irritation to tongue and oral mucosa.

Profile Diameter	Ref.C
Ø4.0	AANOHC4000T
Ø5.0	IHC400T
Ø6.0	AANOHC6000T



Impression Coping (Transfer)


- Guide Pin(AAOTGP10 / AAOTGP12) included
- Should be tightened with Impression Coping Driver (Page.078)
- Special impression coping screw which can be used with a 1.2mm hex driver is available on request.

Profile Diameter	Height (mm)	Type	Ref.C
Ø4.0	7.5	Octa	AAOITO4010T
		Non-octa	AAOITN4010T
	9.5	Octa	AAOITO4012T
		Non-octa	AAOITN4012T
Ø5.0	7.5	Octa	AAOITO5010T
		Non-octa	AAOITN5010T
	9.5	Octa	AAOITO5012T
		Non-octa	AAOITN5012T
Ø6.0	7.5	Octa	AAOITO6010T
		Non-octa	AAOITN6010T
	9.5	Octa	AAOITO6012T
		Non-octa	AAOITN6012T

Impression Coping (Pick-up)


- Guide Pin included

Profile Diameter	Height (mm)	Type	Ref.C
Ø4.0	10.0	Octa	AAOIPO4010T
		Non-octa	AAOIPN4010T
	12.0	Octa	AAOIPO4012T
		Non-octa	AAOIPN4012T
Ø5.0	10.0	Octa	AAOIPO5010T
		Non-octa	AAOIPN5010T
	12.0	Octa	AAOIPO5012T
		Non-octa	AAOIPN5012T
Ø6.0	10.0	Octa	AAOIPO6010T
		Non-octa	AAOIPN6010T
	12.0	Octa	AAOIPO6012T
		Non-octa	AAOIPN6012T

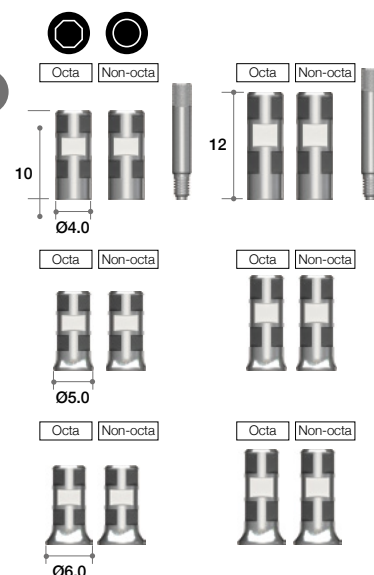

Lab Analog

Profile Diameter	Ref.C
Ø3.8	AANOLA4000
Ø4.8	IOA300
Ø5.8	AANOLA6000

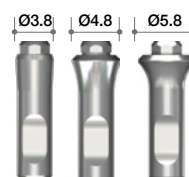


Temporary Cylinder

- Cylinder Screw (IRCS200) included

• Recommend torque : 25Ncm

Profile Diameter	Type	Ref.C
Ø4.0	Octa	AANOTCO4010T
	Non-octa	AANOTCN4010T
Ø5.0	Octa	AANOTCO5010T
	Non-octa	AANOTCN5010T
Ø6.0	Octa	AANOTCO6010T
	Non-octa	AANOTCN6010T

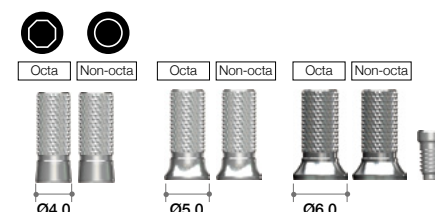


EZ Post Cylinder

- Cylinder Screw (IRCS200) included

• Recommend torque : 25Ncm

Profile Diameter	Post Height(mm)	Type	Ref.C
Ø4.0	5.5	Octa	AAOECO4005T
	7.0		AAOECO4007T
	5.5	Non-octa	AAOECN4005T
	7.0		AAOECN4007T
Ø5.0	5.5	Octa	AAOECO5005T
	7.0		AAOECO5007T
	5.5	Non-octa	AAOECN5005T
	7.0		AAOECN5007T
Ø6.0	5.5	Octa	AAOECO6005T
	7.0		AAOECO6007T
	5.5	Non-octa	AAOECN6005T
	7.0		AAOECN6007T

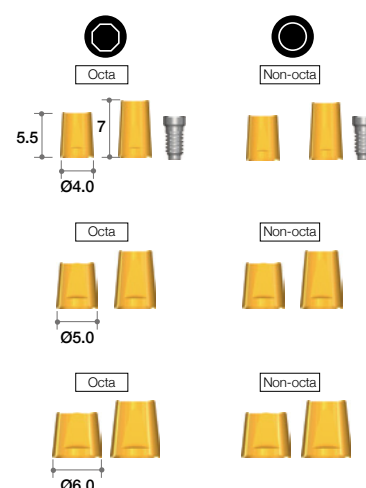


➔ Components for Octa Abutment

Gold Cylinder

- Cylinder Screw (IRCS200) included

- For customizing abutment for screw retained multi-unit restoration.
- Available in both octa(red) and non-octa(white)
- Melting point of gold alloy : 1063°C
- Threaded sleeves allow for better retention of resin or wax.
- Available in three diameters (Ø4.0, Ø5.0 & Ø6.0).
- Recommend torque : 30Ncm

Profile Diameter	Type	Ref.C
Ø4.0	Octa	AANGCO4000T
	Non-octa	AANGCN4000T
Ø5.0	Octa	ILOGO100T
	Non-octa	ILOGN100T
Ø6.0	Octa	AANGCO6000T
	Non-octa	AANGCN6000T

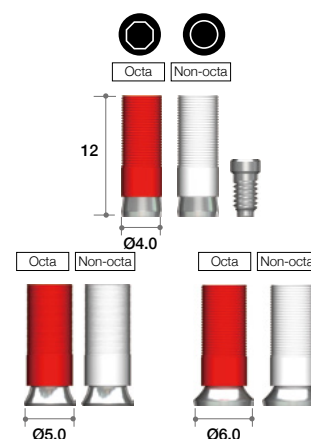


CCM Cylinder

- Cylinder Screw (IRCS200) included

- Useful to make a customized abutment in difficult situations.
- Can be casted with non-precious alloys (Ni-Cr, Cr-Co alloys).
- Non-precious melting temperature : Depend on Manufacturer.
- Threaded sleeves for convenient Resin/ Wax-up.
- Melting temperature of CCM : 1300~1400°C
- Recommend torque : 35Ncm

Profile Diameter	Type	Ref.C
Ø4.0	Octa	AANCCO4000T
	Non-octa	AANCCN4000T
Ø5.0	Octa	AANCCO5000T
	Non-octa	AANCCN5000T
Ø6.0	Octa	AANCCO6000T
	Non-octa	AANCCN6000T

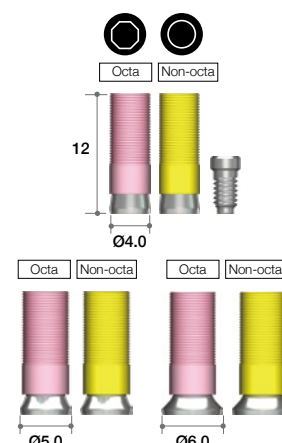


Plastic Cylinder

- Cylinder Screw (IRCS200) included

- Economical option
- Used for customizing abutment for screw retained multi-unit restorations.
- Available in both octa(red) and non-octa(white)
- Threaded sleeves allow for better retention of resin or wax.
- Recommend torque : 25Ncm

Profile Diameter	Type	Ref.C
Ø4.0	Octa	AAOTCO4010T
	Non-octa	AAOTCN4010T
Ø5.0	Octa	IOPH100T
	Non-octa	IOPN100T
Ø6.0	Octa	AAOTCO6010T
	Non-octa	AAOTCN6010T

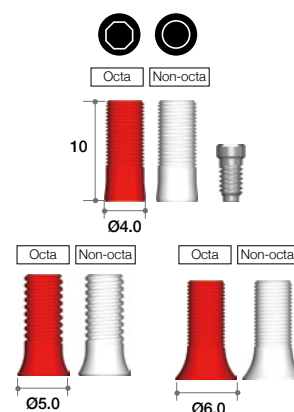


➡ Components for Octa Abutment_Digital (Continued)

Scan Abutments

- Abutment Screw(SAIRCS200) included.

- For Chairside/ Labside
- Supporting Dental CAD
 - 3Shape / Exocad / Dental Wings
- Recommend torque : By Hand (5~8Ncm)

Profile Diameter	Height (mm)	Ref.C
Ø4.0	11	AOCES4011T

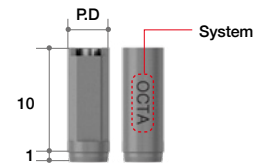


Analog

- Analog Screw(ALS18) included.

- For Chairside/ Labside
- Supporting Dental CAD
 - 3Shape
 - exocad
- 2 piece type


Profile Diameter	Ref.C
Ø3.8	OCTAALST
Ø4.8	OCTAALRT
Ø5.8	OCTAALWT

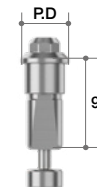

ZrGEN Abutments

- Abutment Screw(IRCS200) included.

- Titanium Base
- 1 set consists of 10 Abutments.
 - included spare Abutment Screw
- Supporting Dental CAD
 - 3Shape
 - Exocad
 - Dental Wing
- Post Height can be checked by the number of Groove.
 - Post Height : 5mm → Groove : 2ea
 - Post Height : 6mm → Groove : 4ea
 - Post Height : 8mm → Groove : 6ea
- Recommend torque : 35Ncm

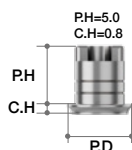

System	Profile Diameter	Cuff Height (mm)	Post Height (mm)	Type	Ref.C
Octa Level	Small	Ø5.0	0.8	Octa	AOCEPS5015.MTN
			6		AOCEPS5016.MTN
			8		AOCEPS5018.MTN
		Ø5.0	0.8	Non-Octa	ANOEPS5015.MTN
			6		ANOEPS5016.MTN
			8		ANOEPS5018.MTN
	Regular	Ø5.5	0.8	Octa	AOCEPR5515.MTN
			6		AOCEPR5516.MTN
			8		AOCEPR5518.MTN
		Ø5.5	0.8	Non-Octa	ANOEPW5515.MTN
			6		ANOEPW5516.MTN
			8		ANOEPW5518.MTN
	Wide	Ø6.5	0.8	Octa	AOCEPW6515.MTN
			6		AOCEPW6516.MTN
			8		AOCEPW6518.MTN
		Ø6.5	0.8	Non-Octa	ANOEPW6515.MTN
			6		ANOEPW6516.MTN
			8		ANOEPW6518.MTN

➔ Components for Octa Abutment_Digital

TiGEN Abutments

(MegaGen type)

- Abutment Screw(IRCS200) included.

- Pre-milled Abutment
- 1 set consists of 10 Abutments.
 - included spare Abutment Screw
- Supporting DentalCAD
 - 3Shape / exocad / Dental Wings
- Supporting Milling Machine
 - MegaGen Implant : BX5
 - ARUM DENTISTRY
- Recommend torque : 35Ncm

System	Color	Profile Diameter	Height (mm)	Type	Ref.C
Octa Level	Small	Silver	20	Octa	OCTS1020.MTN
				Non-Octa	NOTS1020.MTN
				Octa	OCTS1220.MTN
				Non-Octa	NOTS1220.MTN
	Octa			OCTR1020.MTN	
	Non-Octa			NOTR1020.MTN	
	Octa			OCTR1220.MTN	
	Non-Octa			NOTR1220.MTN	
	Regular	Ø10		Octa	OCTW1020.MTN
				Non-Octa	NOTW1020.MTN
				Octa	OCTW1220.MTN
				Non-Octa	NOTW1220.MTN
Wide	Ø12	Octa	OCTW1020.MTN		
		Non-Octa	NOTW1020.MTN		
		Octa	OCTW1220.MTN		
		Non-Octa	NOTW1220.MTN		




(NT type)

- Abutment Screw(IRCS200) included.

- Pre-milled Abutment
- 1 set consists of 10 Abutments
 - included spare Abutment Screw
- Supporting DentalCAD
 - 3Shape / exocad
- Recommend torque : 35Ncm

System	Color	Profile Diameter	Height (mm)	Type	Ref.C
Octa Level	Small	Silver	16	Octa	OCTSN1016.MTN
				Non-Octa	NOTSN1016.MTN
				Octa	OCTSN1216.MTN
				Non-Octa	NOTSN1216.MTN
	Octa			OCTRN1016.MTN	
	Non-Octa			NOTRN1016.MTN	
	Octa			OCTRN1216.MTN	
	Non-Octa			NOTRN1216.MTN	
	Regular	Ø10		Octa	OCTWN1016.MTN
				Non-Octa	NOTWN1016.MTN
				Octa	OCTWN1216.MTN
				Non-Octa	NOTWN1216.MTN
Wide	Ø12	Octa	OCTWN1016.MTN		
		Non-Octa	NOTWN1016.MTN		
		Octa	OCTWN1216.MTN		
		Non-Octa	NOTWN1216.MTN		

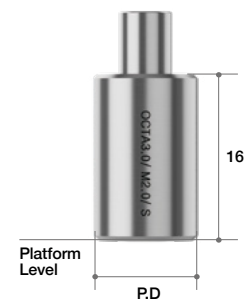

(Medentika type)

- Abutment Screw(IRCS200) included.


- Pre-milled Abutment
- 1 set consists of 10 Abutments
 - included spare Abutment Screw
- Supporting DentalCAD
 - 3Shape / exocad
- Recommend torque : 35Ncm

System		Color	Profile Diameter	Height (mm)	Type	Ref.C
Octa Level	Small	Silver	Ø12	14	Octa	OCTSM1214.MTN
					Non-Octa	NOTSM1214.MTN
	Regular				Octa	OCTRM1214.MTN
					Non-Octa	NOTRM1214.MTN
	Wide				Octa	OCTWM1214.MTN
					Non-Octa	NOTWM1214.MTN



II. Abutment Level Prosthesis

3-1. Multi-unit Abutment & Components (All-on-4) (N-Type)

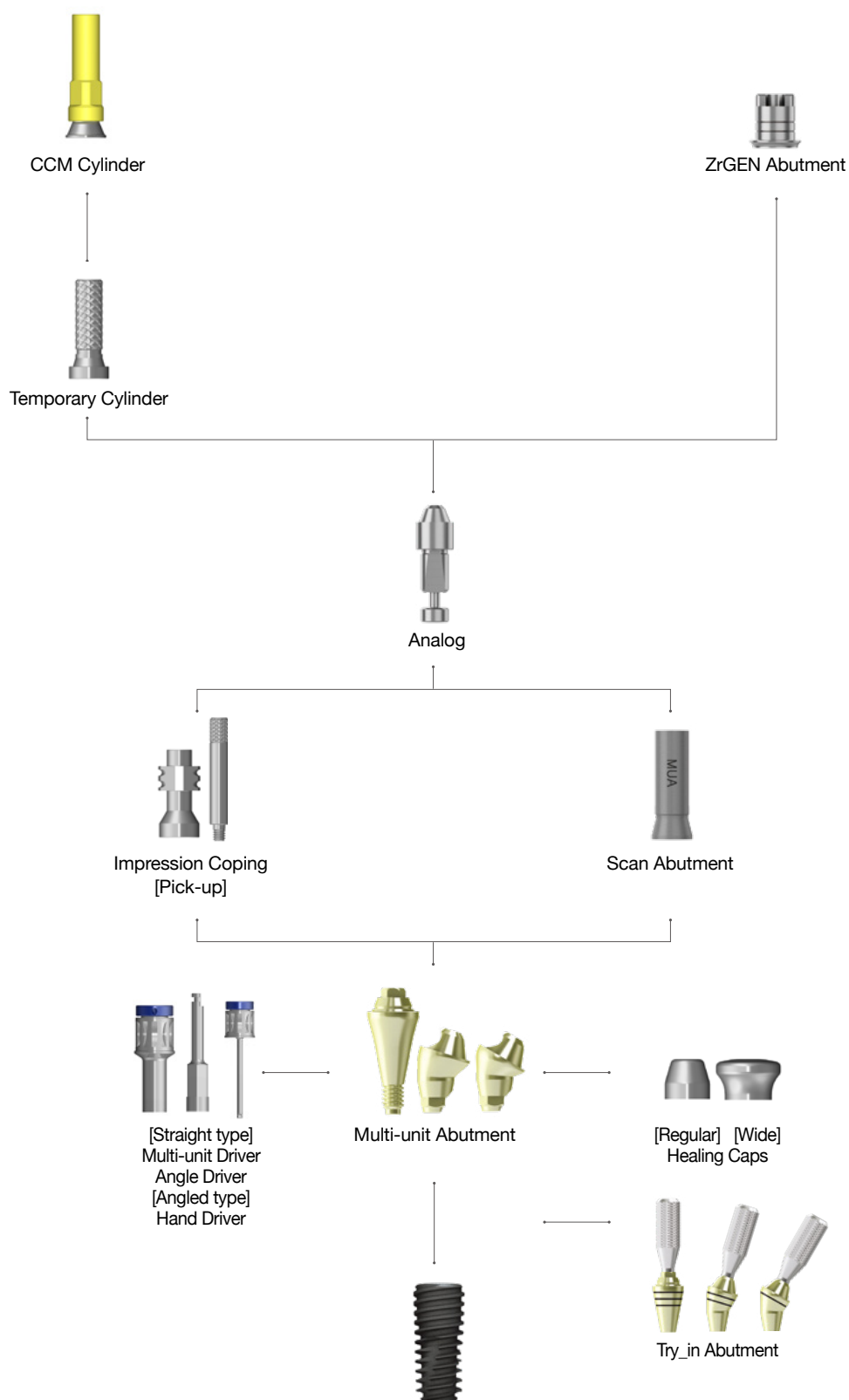


➔ Multi-unit Abutment

Multi-unit Abutment - Straight

- MUA Straight Carrier (MUASC) included

- Regular: Compatible with AnyOne Internal System
- Recommend torque : Mini / 25Ncm, Regular / 35Ncm
- Mini Connection : Coming Soon

Cuff Height (mm)	Connection	Type	Ref.C
1.5	Mini	1-piece (M2)	MUASTMN0015C
2.5			MUASTMN0025C
3.5			MUASTMN0035C
4.5			MUASTMN0045C
1.5	Regular		MUAAON0015C
2.5			MUAAON0025C
3.5			MUAAON0035C
4.5			MUAAON0045C

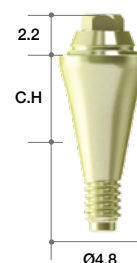


Multi-unit Angled Abutment - 17°

- MUA Screw (MUAOS) included
- MUA Angled Carrier (MUAAC) included

- Recommend torque : Mini / 25Ncm, Regular / 25Ncm
- Mini Connection : Coming Soon

Cuff Height (Labial) (mm)	Connection	Ref.C
2.5 (1.1)	Mini	MUASTMH1725TC
3.5 (2.1)		MUASTMH1735TC
4.5 (3.1)		MUASTMH1745TC
2.5 (1.1)		MUASTMN1725TC
3.5 (2.1)		MUASTMN1735TC
4.5 (3.1)		MUASTMN1745TC
2.5 (1.1)	Regular	MUAAOH1725TC
3.5 (2.1)		MUAAOH1735TC
4.5 (3.1)		MUAAOH1745TC
2.5 (1.1)		MUAAON1725TC
3.5 (2.1)		MUAAON1735TC
4.5 (3.1)		MUAAON1745TC

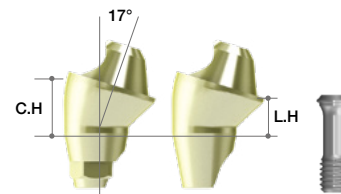


Multi-unit Angled Abutment - 30°

- MUA Screw (MUAOS) included
- MUA Angled Carrier (MUAAC) included

- Recommend torque : Mini / 25Ncm, Regular / 25Ncm
- Mini Connection : Coming Soon

Cuff Height (Labial) (mm)	Connection	Ref.C
3.5 (1.1)	Mini	MUASTMH3035TC
4.5 (2.1)		MUASTMH3045TC
3.5 (1.1)		MUASTMN3035TC
4.5 (2.1)		MUASTMN3045TC
3.5 (1.1)	Regular	MUAAOH3035TC
4.5 (2.1)		MUAAOH3045TC
3.5 (1.1)		MUAAON3035TC
4.5 (2.1)		MUAAON3045TC

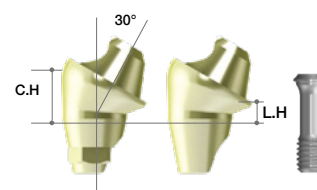

►► Contents of Multi-unit Abutment Set

Multi-unit Abutment Healing Cap-type Set Reference Code


Order code : Available by changing to 'HP' instead of 'C' or 'TC' from current Ref.C

Ex) MUAAOH1725TC → MUAAOH1725 **HP**

Multi-unit Abutment CCM-type Set Reference Code

Order code : Available by changing to 'P' instead of 'C' or 'TC' from current Ref.C

Ex) MUAAOH1725TC → MUAAOH1725 **P**

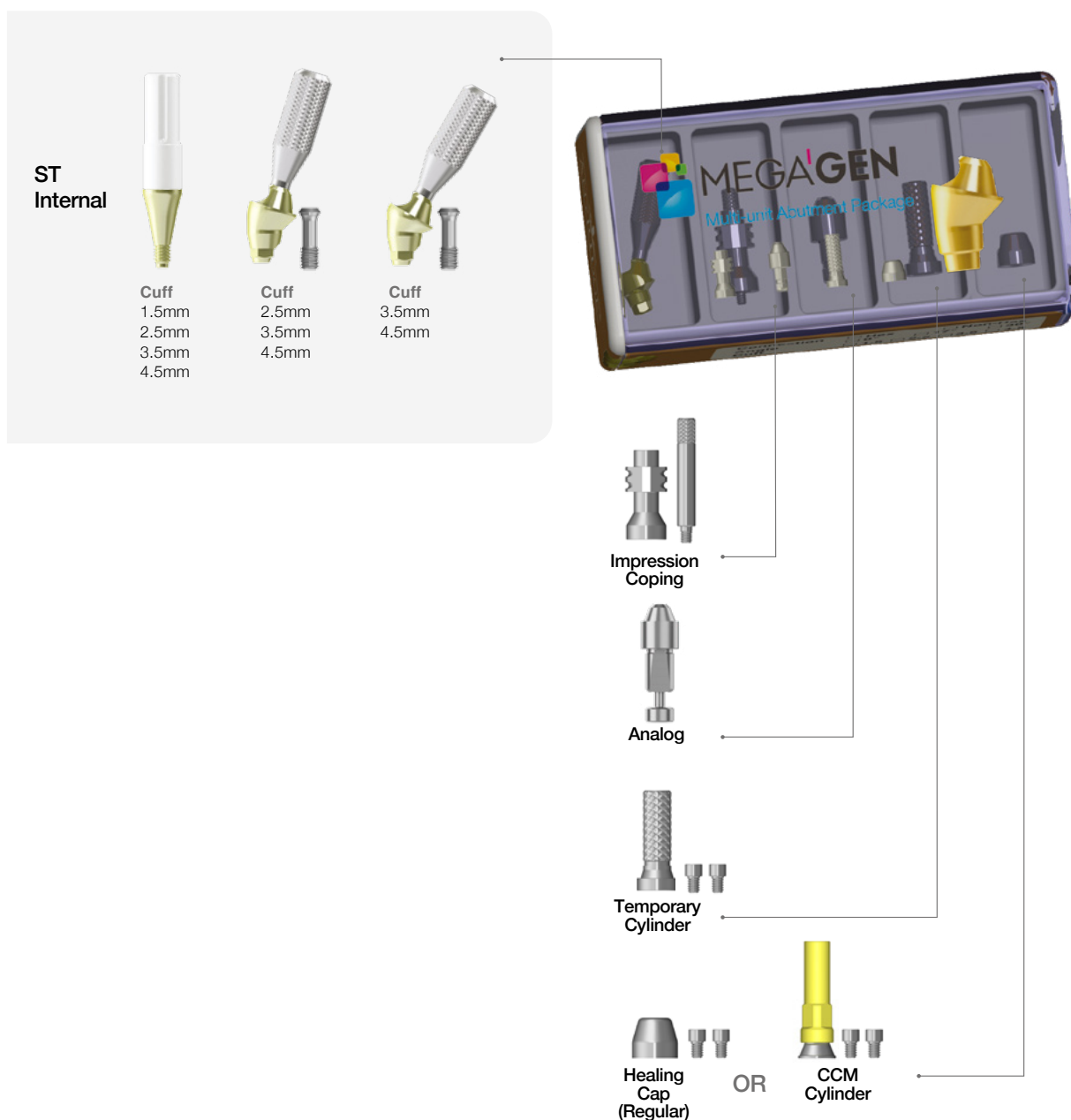


➡ Components for Multi-unit Abutment (Continued)

Try-in Abutment

- Cuff height is indicated with laser markings
- Straight, 17°, 30°
- Mini Connection : Coming Soon

Angle	Connection	Cuff Marking	Ref.C
Straight	Mini	1.5 / 2.5 / 3.5 / 4.5	MUTIASTM00C
17°		2.5 / 3.5 / 4.5	MUTIASTM17C
30°		3.5 / 4.5	MUTIASTM30C
Straight	Regular	1.5 / 2.5 / 3.5 / 4.5	MUTIAAO00C
17°		2.5 / 3.5 / 4.5	MUTIAAO17C
30°		3.5 / 4.5	MUTIAAO30C

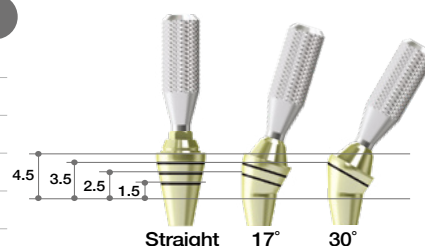

Try-in Abutment Set reference code

Order code : MUTIAAO00C P


* Kit contains Straight, 17° and 30° type of Try-in Abutments (1 each)



Healing Cap

- Cylinder Screw (MUAS) 2ea included
- The size of healing cap can be selected depending on soft tissue volume or type of restorations.

Type	Ref.C
Regular	MUAHCL
Wide	MUAHCWL

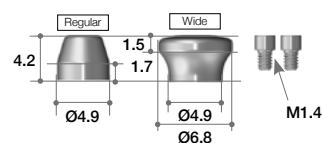

Healing Cap Set reference code

Order code : Available by changing to 'P' instead of 'L' from current Ref.C



Ex) MUAHCL → MUAHCP



▶▶ Starting Package Contents



	Type	Ref.C
Healing Cap	Hex	SKAOHN3000H
	Non Hex	SKAONN3000H
CCM Abutment	Hex	SKAOHN3000
	Non Hex	SKAONN3000

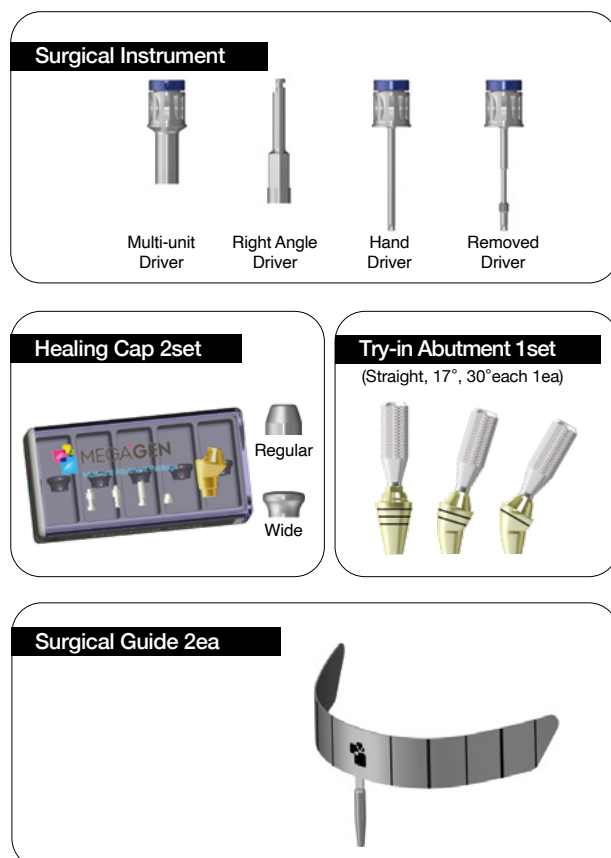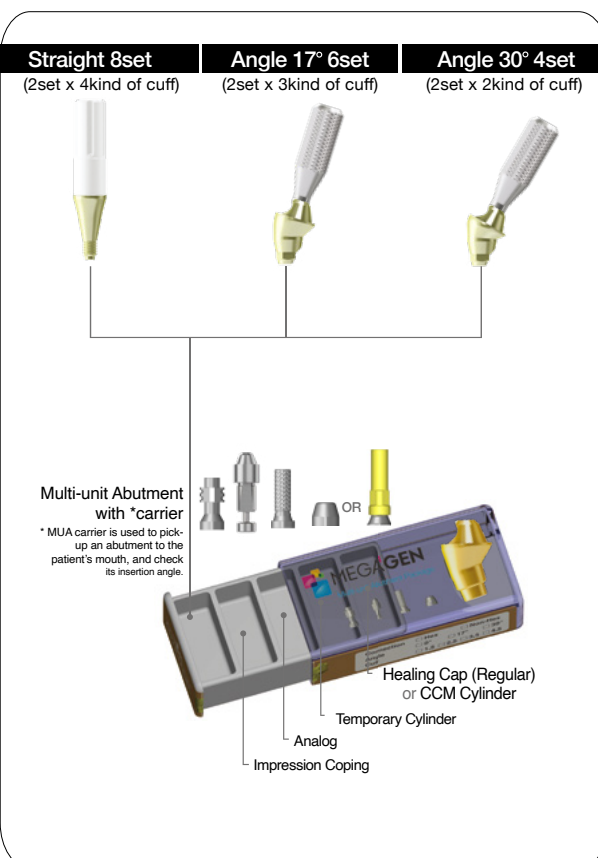


➡ Components for Multi-unit Abutment (Continued)

Impression Coping

(Pick-up)

- Guide pin (MUAGP) included
- For use with impression taking at abutment level


Connection	Ref.C
Non-Hex	MUAICT

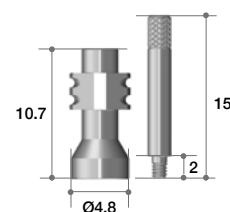

Analog

- For use with duplicating multi-unit abutment in working model
- Available as RP Analog for 3D-printed working model


Head form	Ref.C
Multi-unit Abutment(Nobel)	MUAALT

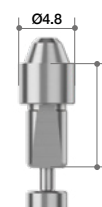

Temporary Cylinder

- Cylinder screws (MUAS) 2EA included
- For use with fabricating acrylic provisional restoration
- Grooves on post cylinder for storing resin adhesion
- Back-up screw is included
- Recommended torque: 15Ncm


Connection	Ref.C
Non-Hex	MUATCL

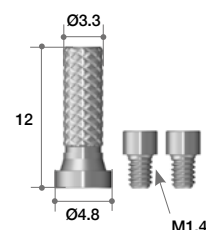

CCM Cylinder

- Cylinder screws (MUAS) 2EA included
- For use with fabricating screw-retained prosthesis with metal-reinforced or bar-structured overdenture
- Can be cast using non-precious dental alloys (Ni-Cr, Cr-Co alloys)
- Melting temperature of CCM base: 1300~1400°C
- Back-up screw included
- Recommended torque: 15Ncm


Connection	Ref.C
Non-Hex	MUACCML

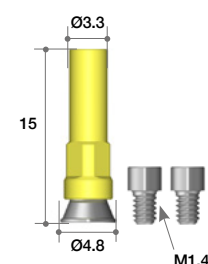

Scan Abutments

- Abutment Screw(SAMUAS) included
- For Chairside/ Labside
- Supporting Dental CAD
 - 3Shape / Exocad / Dental Wings
- Recommend torque : By Hand (5~8Ncm)

Profile Diameter	Height (mm)	Ref.C
Ø4.0	13	AMUASR4013T

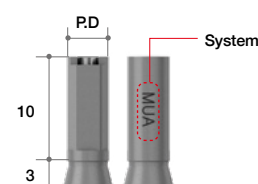


ZrGEN Abutments

- Cylinder Screw(MUAS) included.

- Titanium Base
- 1 set consists of 10 Abutments.
 - included spare Abutment Screw
- Supporting Dental CAD
 - 3Shape / Exocad / Dental Wing
- Post Height can be checked by the number of Groove.
 - Post Height : 5mm → Groove : 2ea
 - Post Height : 6mm → Groove : 4ea
 - Post Height : 8mm → Groove : 6ea
- Recommend torque : 15Ncm


Profile Diameter	Cuff Height (mm)	Post Height (mm)	Type	Ref.C
Ø5.5	0.8	5	N Type	AMUAPR5515N.MTN
		6		AMUAPR5516N.MTN
		8		AMUAPR5518N.MTN

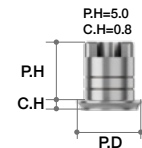

Multi-unit Driver

- Use to torque straight type Multi-unit Abutments.
- Use with a torque wrench (**ref code: MTW300A**)


Hex	Length	Ref.C
2.0	10	MUD10

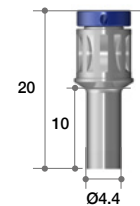

Right Angle Driver

- Use to torque straight type Multi-unit Abutments.
- Use with latch-type handpiece.
- Use with Meg-TORQ (**ref code: MEG_TORQ**)


Hex	Length	Ref.C
2.0	10	MURAD10

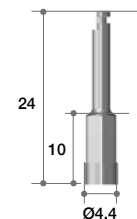

Hand Driver

- Use for abutment screw with 1.2 hex hole.
- Use up to 15° divergent.
- It should use under 30Ncm torque.

Hex	Length	Ref.C
1.2	20	MUHD1220

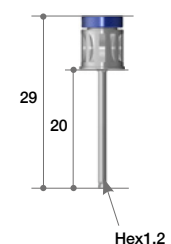


➡ Components for Multi-unit Abutment

Removal Driver

- Use for abutment screw with 1.2 hex hole.
- Use up to 15° divergent.
- Exclusively for AnyRidge system.
- It should use under 30Ncm torque.


Hex	Length	Ref.C
1.2	20	MUARD20

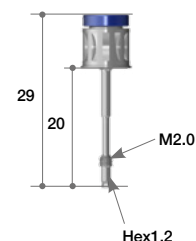

Surgical Guide


- 7mm distance between lines
- Place center pin after initial drilling at center of arch

Angle	Marking Length (mm)	Ref.C
30°	7	MUSG70



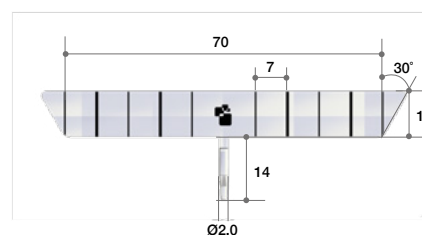


►► How to use Surgical Guide

- * For easy identification, surgical guide includes thicker lines for canines & second premolars, as most common indicators
- * Surgical guide can also be used with first molars

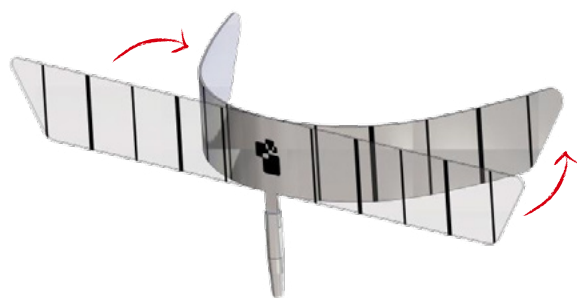


Bend to use

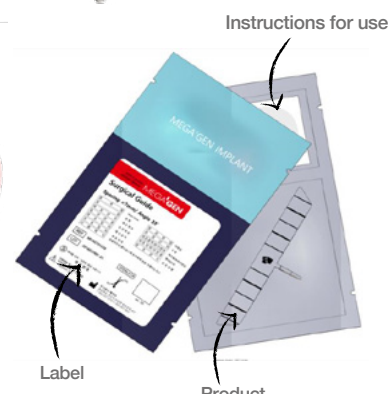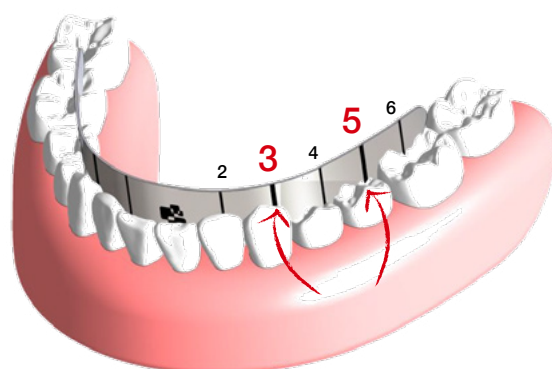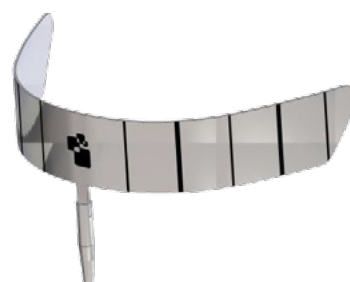

[Packaging]

II. Abutment Level Prosthesis

3-2. Multi-unit Abutment & Components (All-on-4) (S-Type)

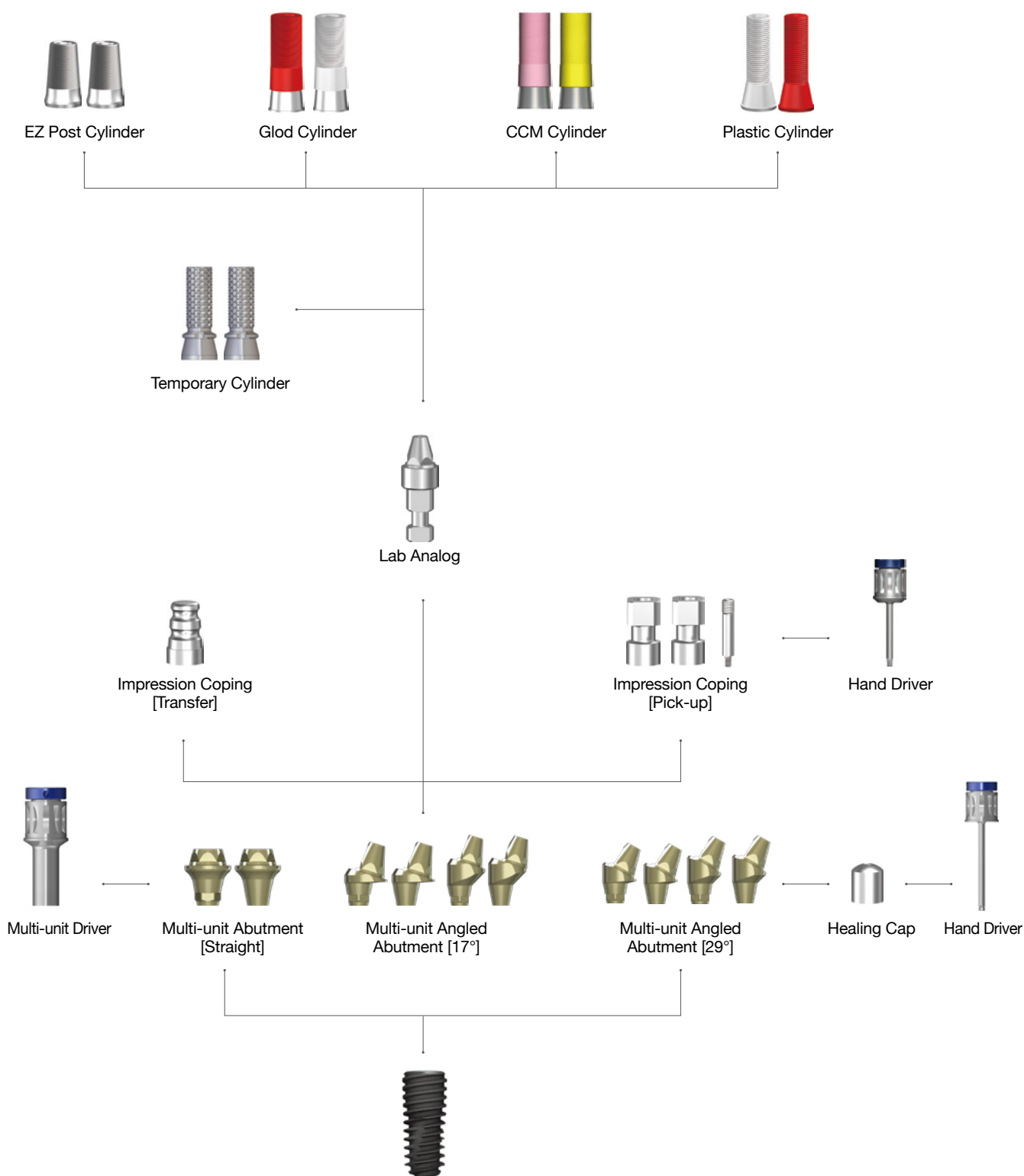


➔ Components for Multi-unit Abutment (Continued)

Multi-unit Abutment (Straight)

- Multi-unit Abutment Screw (MUS15 / MUS25 / MUS35 / MUS45 / MUS55) included.
- Recommend torque : 35Ncm


Cuff Height (mm)	Connection	Ref.C
1.5	Regular	Hex
2.5		
3.5		
4.5		
5.5		
1.5	Regular	Non-Hex
2.5		
3.5		
4.5		
5.5		

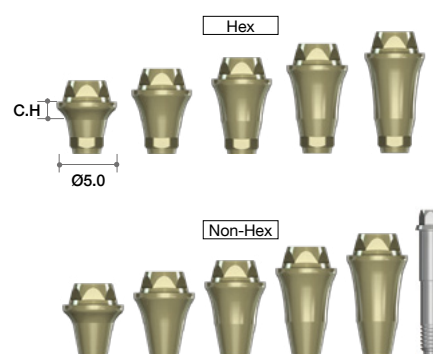

Multi-unit Angled Abutment (17°)

- Abutment Screw (MUAS20) included
- Recommend torque : 35Ncm


Cuff Height (mm)	Connection	Ref.C
1.0	Regular	Hex
2.0		
3.0		
4.0		
1.0	Regular	Non-Hex
2.0		
3.0		
4.0		

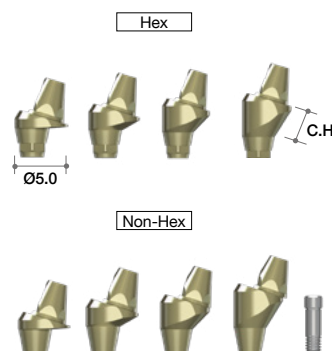

Multi-unit Angled Abutment (29°)

- Abutment Screw (MUAS20) included
- Recommend torque : 35Ncm

Cuff Height (mm)	Connection	Ref.C
1.0	Regular	Hex
2.0		
3.0		
4.0		
1.0	Regular	Non-Hex
2.0		
3.0		
4.0		

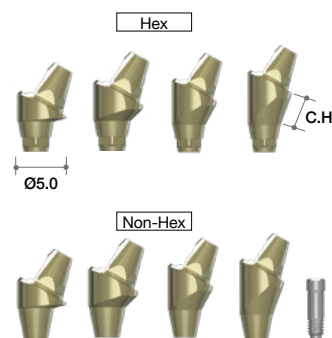

Healing Cap

Profile Diameter	Ref.C
Ø5.0	REC600



Impression Coping (Transfer)

Profile Diameter	Ref.C
Ø4.8	RITE480



Impression Coping (Pick-up)

- Guide Pin (RICG150) included

Height (mm)	Ref.C
9.4	RIEH480T
	RIEN480T



Lab Analog


Profile Diameter	Ref.C
Ø4.8	RELA300



Temporary Cylinder

- Cylinder Screw (TASH140) included
• Recommend torque : 15Ncm

Profile Diameter	Ref.C
Ø4.8	ETH100T
	ETN100T

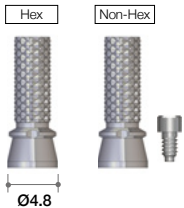


➔ Components for Multi-unit Abutment

EZ Post Cylinder

- Cylinder Screw (TASH140) included
- Recommend torque : 15Ncm

Profile Diameter	Type	Ref.C
Ø5.0	Hex	RCA900T
	Non-Hex	RCA800T

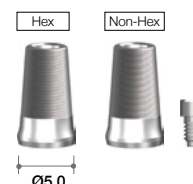

Gold Cylinder

- Cylinder Screw (TASH140) included
- For customizing abutment for screw retained multi-unit restoration.
 - Available in both octa(red) and non-octa(white)
- Melting point of gold alloy : 1063°C
- Threaded sleeves allow for better retention of resin or wax.
- Recommend torque : 15Ncm


Profile Diameter	Sleeve color	Ref.C
Ø4.8	Red	REGC200T
	White	REGC100T



CCM Cylinder

- Cylinder Screw (TASH140) included
- Useful to make a customized abutment in difficult situations.
- Can be casted with non-precious alloys (Ni-Cr, Cr-Co alloys).
- Non-precious melting temperature : Depends on Manufacturer
- Threaded sleeves for convenient Resin/ Wax-up.
- Melting temperature of CCM : 1300~1400°C
- Recommend torque : 15Ncm

Profile Diameter	Sleeve color	Ref.C
Ø4.8	Pink	RCA5013HT
	Yellow	RCA5013NT

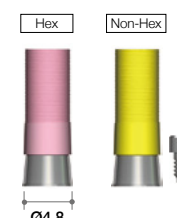


Plastic Cylinder

- Cylinder Screw (TASH140) included
- Economical option
- Used for customizing abutment for screw retained multi-unit restorations.
 - Available in both Hex(red) and Non-Hex(white)
- Threaded sleeves allow for better retention of resin or wax.
- Recommend torque : 15Ncm


Profile Diameter	Sleeve color	Ref.C
Ø5.2	Red	RPEH100T
	White	RPEN100T



- For cutting & shaping cortical bone to secure initial stability in hard bone (D1 & D2)

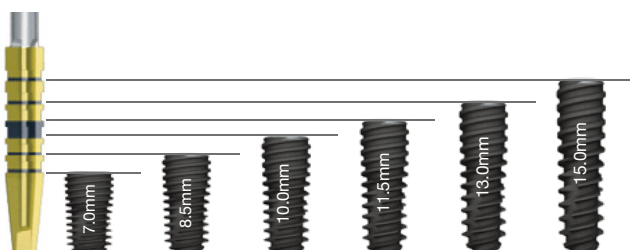


► ► Surgical Drilling Sequence

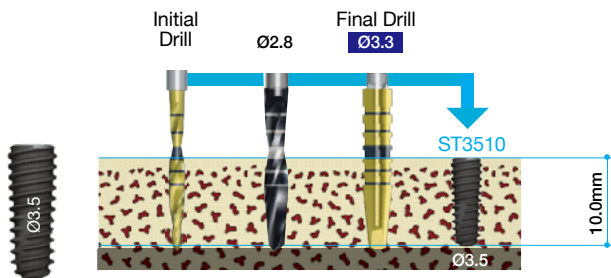

Ø3.5 Fixture

Ø3.5 Drilling Sequence


Actual fixture length is 10.0mm & shaping drill length is 0.59mm, so total drill depth is 10.59mm

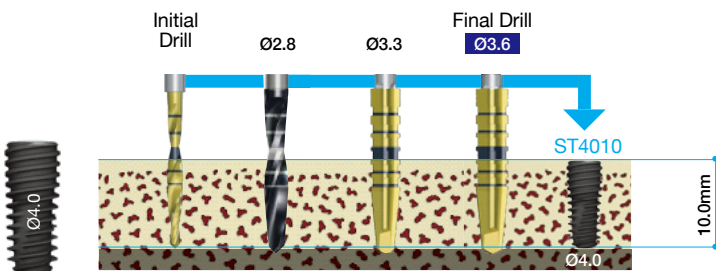

Ø4.0 Fixture

Ø4.0 Drilling Sequence


Actual fixture length is 10.0mm & shaping drill length is 0.68mm, so total drill depth is 10.68mm

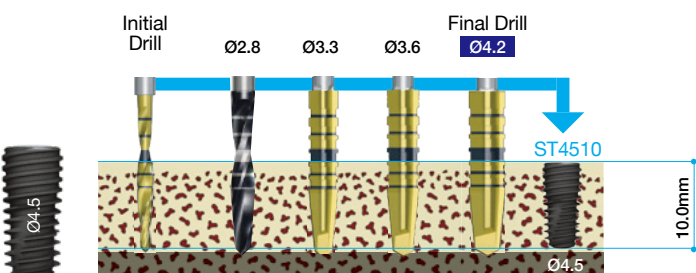

Ø4.5 Fixture

Ø4.5 Drilling Sequence

Actual fixture length is 10.0mm & shaping drill length is 0.85mm, so total drill depth is 10.85mm

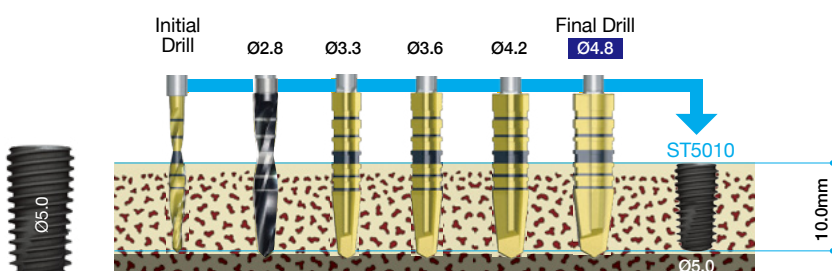

Ø5.0 Fixture

Ø5.0 Drilling Sequence

Actual fixture length is 10.0mm & shaping drill length is 0.89mm, so total drill depth is 10.89mm


➡ ST Surgical Kit Components (Continued)

Lance Drill

- For easy marking of cortical bone for accurate drilling

Diameter	Type	Ref.C
Ø2.0	Long	MGD100L
	Short	*LD2015
	Long	*LD2025
	Ultra-Long	*LD2030


(*) Separate sales item.

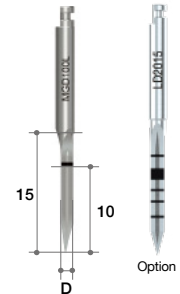

Lindermann Drill

- For controlling drilling path

Diameter	Ref.C
Ø2	TEEL200M

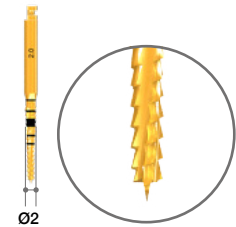

Shaping Drill

- Depth markings up to 7.0~15.0mm
- Dual marking system (groove & laser markings) designed for easy recognition of placement depth
- Coated with TiN to resist corrosion & damage

Diameter	Length(mm)	Ref.C
Ø2.0	33	SD2018S
	38	*SD2018M
	43	*SD2018L
Ø2.8	33	SD2818S
	38	*SD2818M
	43	*SD2818L
Ø3.3	33	SD3318S
	38	*SD3318M
	43	*SD3318L
Ø3.6	33	SD3618S
	38	*SD3618M
	43	*SD3618L
Ø4.2	33	SD4218S
	38	*SD4218M
	43	*SD4218L
Ø4.8	33	SD4818S
	38	*SD4818M
	43	*SD4818L

(*) Separate sales item.



Dense Drill

- For cutting & shaping cortical bone to secure initial stability in hard bone (D1, D2)
- Coated with TiN to resist corrosion & damage

Diameter	Type	Ref.C
Ø3.9	Long	STDD39
Ø4.3		STDD43
Ø4.8		STDD48
Ø5.3		STDD53



➔ ST Surgical Kit Components

Stopper Drill


Diameter	Blade Length (mm)	Ref.C
Ø2.0	7	SD2007M
	8.5	SD2008M
	10	SD2010M
	11.5	SD2011M
Ø2.8	7	SD2807M
	8.5	SD2808M
	10	SD2810M
	11.5	SD2811M
Ø3.3	7	SD3307M
	8.5	SD3308M
	10	SD3310M
	11.5	SD3311M
Ø3.6	7	SD3607M
	8.5	SD3608M
	10	SD3610M
	11.5	SD3611M
Ø4.2	7	SD4207M
	8.5	SD4208M
	10	SD4210M
	11.5	SD4211M
Ø4.8	7	SD4807M
	8.5	SD4808M
	10	SD4810M
	11.5	SD4811M



Flattening Drill

- For flattening irregular bone so that stopper drill can drill to accurate depth
- To use flattening lance, connect to housing (Ø5.0 & Ø6.0) according to final drill diameter
- Ø5.0 = Stopper Drill Ø2.0~ Ø4.3
- Ø6.0 = Stopper Drill Ø4.8~ Ø5.4
- Boundary is formed by housing & guides position of subsequent fixture drilling

Diameter	Length(mm)	Ref.C
Ø5.0 / Ø2.0	3.5	FD5020
Ø6.0 / Ø2.0		FD6020

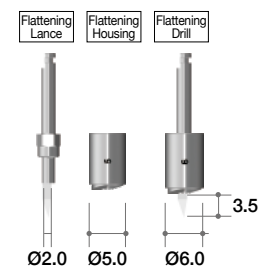

Drill Extension

- For extending length of drill connected to handpiece
- Can be damaged if 45Nm or higher torque is applied


Ref.C
MDE150



Direction Indicator

- For checking location, direction, & depth of drilled hole

Diameter	Ref.C
Ø2.0 / Ø2.8	MDI100
Ø3.0 / Ø4.6	MDI3348

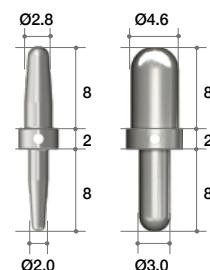

Ratchet Wrench

- For applying stronger force than handpiece
- No bearing system: no blood-triggered corrosion, not easily damaged
- Connect to ratchet extension to use
- Arrow on head indicates direction of force to be applied

Ref.C

MRW040S




Hand Driver (1.2 Hex)

- For all cover screws, abutment screws, & healing abutments
- Various lengths available
- Can be used as torque driver after connecting it to torque wrench without separate adapter
- Hex tip is designed to withstand 35~45Ncm torque

Length(mm)	Type	Ref.C
5	Ultra-short	*TCMHDU1200
10	Short	TCMHDS1200
15	Long	TCMHDL1200
20	Extra-long	*TCMHDE1200
30	Extra-long	*TCMHDE1230

(*) Separate sales item

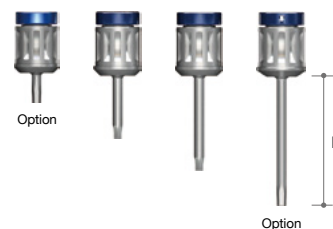

Handpiece Connector

- For retrieving fixture from ampule & positioning in mouth using handpiece
- Spring-type connection allows easier & safer fixture pick-up & positioning
- First (bottom) mark on shaft indicates position of fixture platform, while second (top) mark indicates distance from fixture platform. Lower part indicates 3.5mm from bone level & upper part indicates 4.5mm as guide for flapless surgery.

Length(mm)	Type	Connection	Ref.C
10	Short	Mini	STHCS21
15	Long	Mini	STHCL21
15	Short	Regular	STHCS25
20	Long	Regular	STHCL25

Length(mm)	Type	Connection	Ref.C
15	Short	Mount	STHCS36M



Ratchet Connector

- For placing or removing fixture using ratchet wrench
- Check if ratchet wrench & ratchet connector are connected properly before use
- Too much torque can damage internal hex
- First (bottom) mark on shaft indicates position of fixture platform, while second (top) mark indicates distance from fixture platform. Lower part indicates 3.5mm from bone level & upper part indicates 4.5mm as guide for flapless surgery.

Length(mm)	Type	Connection	Ref.C
10	Short	Mini	STRCS21
15	Long	Mini	STRCL21
15	Short	Regular	STRCS25
20	Long	Regular	STRCL25

Length(mm)	Type	Connection	Ref.C
15	Short	Mount	STRCS36M



Clinical Cases

➡ Clinical Case 1

- Courtesy of Dr. Dae Hee Lee

Fig 1. Due to severe chronic adult periodontitis, bone is quite significantly resorbed across entire upper dentition.

Fig 2. All teeth from premolars to molars were extracted, and waited for more than 2 months.

Fig 3. Lateral sinus graft and GBR were performed on right side.

Fig 4. Same procedures were performed on left side, and perforated sinus membrane was repaired with collagen membrane.

Fig 5. After 5 months, ST implants were placed and bone core biopsies taken from #15 and #25.





Fig 6. Vital bone formation area at 5 mm depth cut > 30%

Fig 7. ST Implant 4010 were placed at #14 and #24, and ST Implant 4510 in other areas.

Fig 8. All ISQ values were good, so impression was taken after 2 months.

Fig 9. Temporary bridge after 2.5 months.

Fig 10. Final prosthesis.

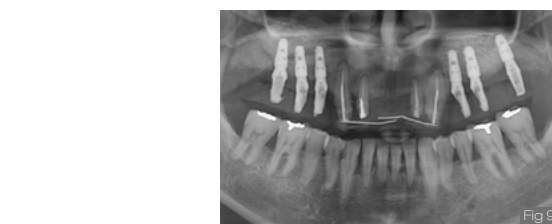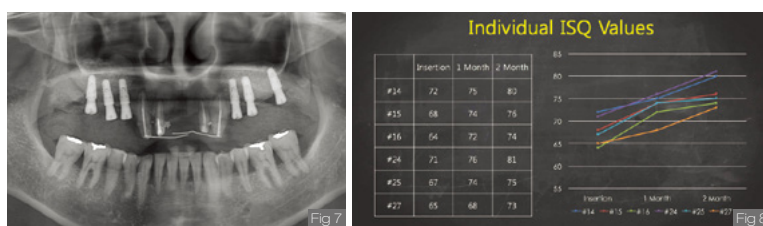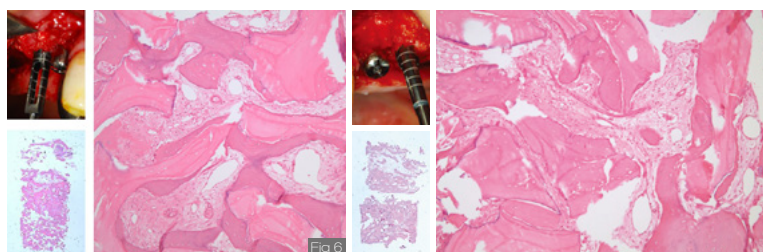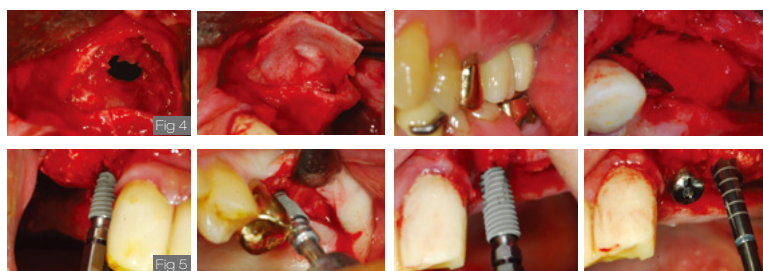

➡ Clinical Case 2


- Courtesy of Dr. Dae Hee Lee

Fig 1. #27 was extracted due to chronic periodontitis, and waited for 2 months.

Fig 2. Two lateral windows were made on septum.

Fig 3. Two ST implants (4510) were placed after 4 months.

Fig 4. 2-year follow-up.

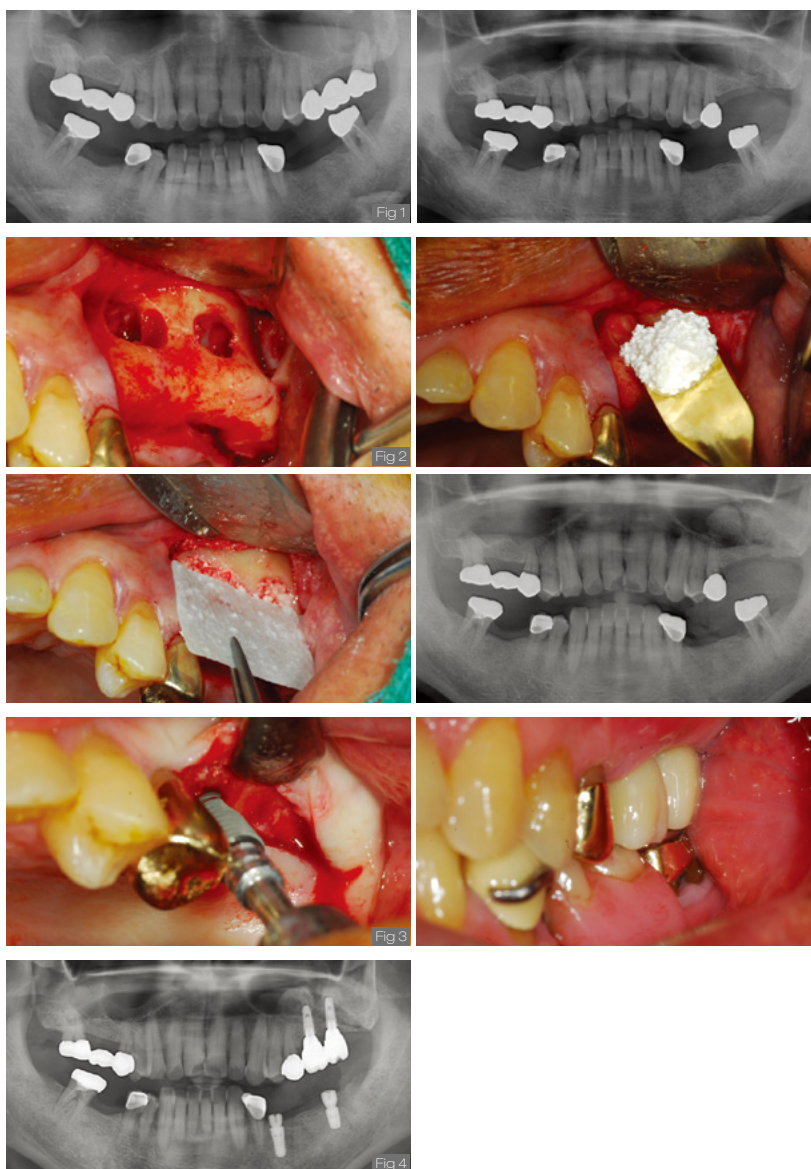



ST™
by MEGA'GEN



București
Timișoara
Iași
Cluj-Napoca

Strada Delea Nouă, Nr. 38, Sector 3
Bvd. Splaiul Nicolae Titulescu, Nr. 8, Parter
Str. Mitropolit Varlaam, Nr. 54, Et. 4
Str. Govora, Nr. 8, Et. 1

Rev. 04