

HU-FRIEDY KEY OPINION LEADER DR. DOMENICO MASSIRONI



MASSIRONI KIT FOR PROSTHETICS BY HU-FRIEDY



Country: **Italy**



Focus on:
Prosthetics



Website:
meg-educational.com



Inspired by:
Prof. Mario Martignoni



Favourite Instrument:
TNMASS1



MASSKITDIN
MASSIRONI KIT FOR
PROSTHETICS



A man who works with his hands is a laborer, a man who works with his hands and his brain is a craftsman; but a man who works with his hands, his brain and his heart is an artist.

Dr. Domenico Massironi




DYNAMIC FUNCTIONALITY IN MODERN PROSTHESIS WITH THE NEW MASSIRONI KIT BY HU-FRIEDY

The new **MASSIRONI KIT** by Hu-Friedy is set of instrumentation used during prosthetic procedures. Conceived by **Dr. Domenico Massironi**, pioneer of the use of the operating microscope in prosthesis, the kit is designed to be used alongside magnification tools, allowing to finalize prosthetic work with maximum precision and accuracy.

The kit consists of 10 Hu-Friedy instruments, including a tissue plier, a high definition mirror as well as a Klemmer plier, which provide a dynamic feature in modern prosthetic dentistry.

The Massironi Kit for prosthetics includes **2 new Enamel Rounded Chisel** instruments, made with **XTS™ technology**, and designed by Dr. Massironi himself in collaboration with Hu-Friedy. These chisels are most often used in the finishing stages of the procedure and the **off-angle design** allows clinicians to access the restoration from any angle.

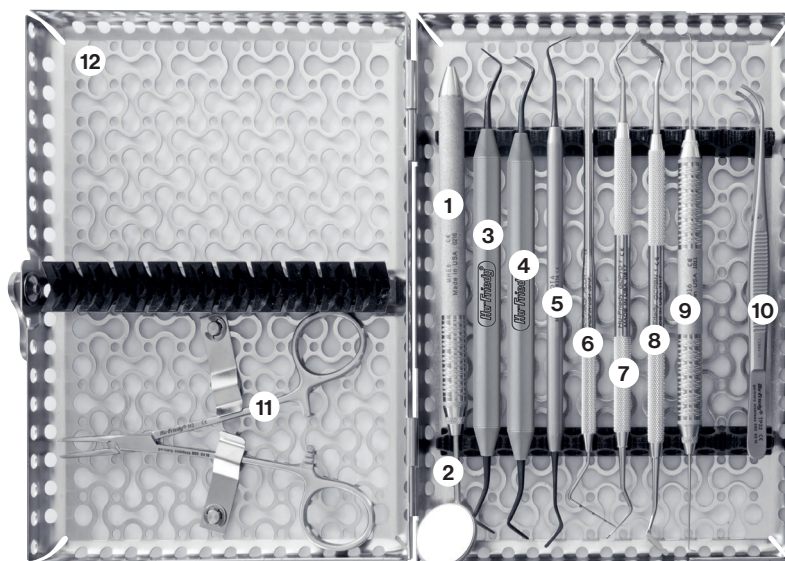

Visit us online at [HuFriedyGroup.eu](https://www.HuFriedyGroup.eu)

MASSIRONI KIT FOR PROSTHETICS BY HU-FRIEDY
HOW TO DO?


THE ENAMEL ROUNDED CHISELS BY MASSIRONI

3 TNMASS1
Enamel Rounded Chisel Massironi, 1,2 mm

The slender width of this chisel makes it ideal for use during highly magnified prosthodontic cases and minimally invasive preparation techniques.

4 TNMASS2
Enamel Rounded Chisel Massironi, 1,5 mm

The slightly larger width of this chisel makes it excellent for everyday use and allows the traditional prosthesis to receive an optimal finishing touch.

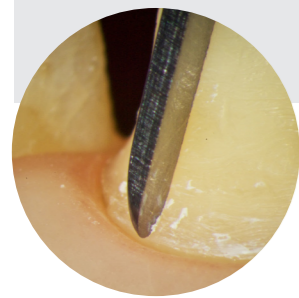

The two chisels are most often used in the finishing stages of the procedure and the off-angle design allows clinicians to access the restoration from any angle.

6 PCORTLINI
Probe Cortellini 0.5-15 hdl #31

Ideal for evaluation of minimum thicknesses in the prosthetic reduction of crowns and laminate veneers, and also for the use of silicone index.




7 GCPYD1
Gingival Cord Packer #1 hdl #41

Spatula for inserting rear deflection wires, also of large diameter, with different angle for optimal work over 360 degrees.



8 GCPBN1
Gingival Cord Packer #BN1 hdl #41

Spatula for inserting front deflection wires, small-medium diameter, with different angle for optimal work over 360 degrees.

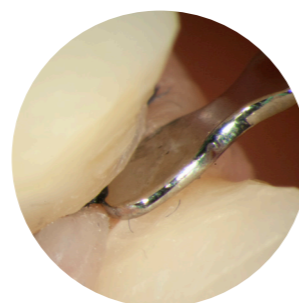

5 CVIPCKOT4
Carver Kotschy #4 interproximal

Multifunctional ultra-thin spatula, suitable for insertion in the deflection phase of thin silk-like threads and for managing the resin in the margin during the relining of the provisional.

9 EXD16
Explorer #1 hdl #6

Ultra-fine speculum of marginal gaps evaluation in try-in and tests of different structures and materials used in prosthetics.

10 TP32
Tissue Pliers Semkin-Taylor #32 anatomical curved 12,5 cm

Small and versatile with perfect grip when using magnifying instruments.

2 MIR5HDE
High Definition Mirror #5 single-sided, European, 1pc./pack

1 MHE6CH
Mirror hdl European hdl #6 Smooth Satin Steel™

Ultra-fine speculum of marginal gaps evaluation in try-in and tests of different structures and materials used in prosthetics.

Used for indirect vision, retraction of lips, cheeks and tongue as well as for illumination of dark areas.

11 H3
Hemostat Halsted-Mosquito #3 curved, 12 cm

Ideal for the removal of provisionals, both front and rear groups.

CREATE YOUR CUSTOM KIT AND CHOOSE THE CASSETTE WITH YOUR FAVORITE COLOR.

Just add the color coding number to the last part of the partcode (replace " _ " with the number)

Red	●	0
Grey	●	1
White	●	2
Orange	●	3
Purple	●	4
Yellow	●	5
Lavender	●	6
Black	●	7
Blue	●	8
Green	●	9
Ocean Blue	●	11

WHY YOU SHOULD USE AN INSTRUMENT MANAGEMENT SYSTEM (IMS):


By standardizing your instruments with IMS, additional time and revenue is generated while extending the useful life of your Hu-Friedy instruments. Dental professionals state that by implementing IMS, they save 5 to 10 minutes per procedure allowing the staff to spend more time on value-added activities.* In addition, safer transportation of soiled instruments reduces risk of injury and exposure to blood borne pathogens.

* When compared to single instrument decontamination. Based on market survey results.

DISCOVER MORE AT IMS WEBPAGE



IME12D1_C2
IMS™ Infinity Series™ Cassette 1/2 DIN for 10 instruments
Dimensions 152 x 184 x 34 mm

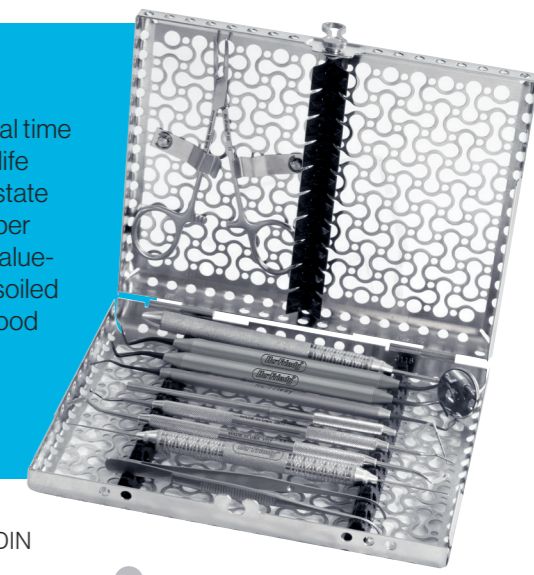

THE MASTER EDUCATIONAL GROUP (MEG)

MELEGNANO - ITALY

TEACHING ACTIVITY:

- 3 Annual courses
- 6 Modular courses
- More than 2,500 hours of lessons
- More than 600 trained dentists

The Master Educational Group (MEG) is an educational center dedicated to esthetic dentistry with a hands-on approach to clinical teaching and education. It offers advanced and continuous trainings with enthusiastic and talented dentists, who share their passion with the participants in exploring methods, innovations and new procedures in a technological operating environment.

Most of the courses and live sessions are held by Dr. Domenico Massironi, who founded the center together with his friend and colleague Dr. Carlo Ghezzi. During his sessions, Dr. Massironi demonstrates his vision of minimally invasive dentistry to maintain tissue integrity as far as possible for harmonious results. The fabrication of a prosthetic restoration consists of many different steps that aim to achieve an optimal anatomical functional and esthetic integration of the prosthetic restoration into the patient's mouth.

www.meg-educational.com/it



For more information or to place an order, call at 00800 48 37 43 39 or visit HuFriedyGroup.eu

Free Call: 00800 48 37 43 39 | Free Fax: 00800 48 37 43 40

E-Mail: info@hu-friedy.eu | Website: HuFriedyGroup.eu

   @HuFriedyGroupEurope

Hu-Friedy Mfg. Co., LLC • European Headquarters • Lyoner Str. 9 • 60528 Frankfurt am Main, Germany • HuFriedyGroup.eu
All company and product names are trademarks of Hu-Friedy Mfg. Co., LLC., its affiliates or related companies, unless otherwise noted.
©2022 Hu-Friedy Mfg. Co., LLC. All rights reserved. HFL-187GB/0722

 **HuFriedyGroup**
The Best In Practice